
















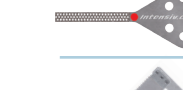


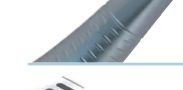





ARTICLE CODES – PRICELIST

	Diamond grit-size	Quantity / Pack	Prices in EUR*
 LAO2018DS/3	Approx Opener DS	Double-sided diamond coated (8µm)	3 pcs 75,80
 LAO2018OS/3	Approx Opener OS	One-sided diamond coated strip (8µm)	3 pcs 75,80
 LIPR-DC-SET	Set IPR-DistanceControl		8 pcs 68,20
 LIPRDC010/3	IPR-DistanceControl	Thickness 0,10 mm	3 pcs 28,00
 LIPRDC015/3	IPR-DistanceControl	Thickness 0,15 mm	3 pcs 28,00
 LIPRDC020/3	IPR-DistanceControl	Thickness 0,20 mm	3 pcs 28,00
 LIPRDC025/3	IPR-DistanceControl	Thickness 0,25 mm	3 pcs 28,00
 LIPRDC030/3	IPR-DistanceControl	Thickness 0,30 mm	3 pcs 28,00
 LIPRDC040/3	IPR-DistanceControl	Thickness 0,40 mm	3 pcs 28,00
 LIPRDC050/3	IPR-DistanceControl	Thickness 0,50 mm	3 pcs 28,00
 LIPRDC100/3	IPR-DistanceControl	Thickness 1,0 mm	3 pcs 28,00
 LOSM1508/6	Ortho-Strips Manuali FE	(8/15µm)	6 pcs 49,45
 LOSM8040/6	Ortho-Strips Manuali CM	(40/80µm)	6 pcs 49,45
 LOS08OP/3	Ortho-Strips Opener DS	Double-sided diamond coated (8µm)	3 pz 76,90
 LOS08OP-L/3	Ortho-Strips Opener OS	One-sided diamond-coated, left (8µm)	3 pz 76,90
 LOS08OP-R/3	Ortho-Strips Opener OS	One-sided diamond-coated, right (8µm)	3 pz 76,90
 LOS15/3	Ortho-Strips	Double-sided diamond coated (15µm)	3 pcs 76,90
 LOS15L/3	Ortho-Strips OS	One-sided diamond-coated, left (15µm)	3 pcs 76,90
 LOS15R/3	Ortho-Strips OS	One-sided diamond-coated, right (15µm)	3 pcs 76,90
 LOS25/3	Ortho-Strips	Double-sided diamond coated (25µm)	3 pcs 76,90
 LOS25L/3	Ortho-Strips OS	One-sided diamond-coated, left (25µm)	3 pcs 76,90
 LOS25R/3	Ortho-Strips OS	One-sided diamond-coated, right (25µm)	3 pcs 76,90
 LOS40/3	Ortho-Strips	Double-sided diamond coated (40µm)	3 pcs 76,90
 LOS40L/3	Ortho-Strips OS	One-sided diamond-coated, left (40µm)	3 pcs 76,90
LOS40R/3	Ortho-Strips OS	One-sided diamond-coated, right (40µm)	3 pcs 76,90
LOS60/3	Ortho-Strips	Double-sided diamond coated (60µm)	3 pcs 76,90
LOS60L/3	Ortho-Strips OS	One-sided diamond-coated, left (60µm)	3 pcs 76,90
LOS60R/3	Ortho-Strips OS	One-sided diamond-coated, right (60µm)	3 pcs 76,90
LOS80C/3	Ortho-Strips Coarse	Double-sided diamond coated (80µm)	3 pcs 76,90
LOS80C-L/3	Ortho-Strips Coarse OS	One-sided diamond-coated, left (80µm)	3 pcs 76,90
LOS80C-R/3	Ortho-Strips Coarse OS	One-sided diamond-coated, right (80µm)	3 pcs 76,90
L071	Ortho-Strips Double-Side	IPR Set refill	6 pcs 153,00
L065L	Ortho-Strips One Sided L	Complete kit of Ortho-Strips diamond on the left side	6 pcs 153,00
L065R	Ortho-Strips One Sided R	Complete kit of Ortho-Strips diamond on the right side	6 pcs 153,00
LOST400	Intensiv Ortho-Strips Tray	Tray delivered with 6 silicone holders	1 pcs 50,80
L056ASS/6	Ricambi Inserti in silicone	1 yellow, 1 orange, 2 white, 1 blue, 1 brown	6 pcs 7,85
LWG-69-A	Intensiv Swingle	Without light	1 pcs 728,00
LWG-69-LT	Intensiv Swingle	With light	1 pcs 1.078,00
L053	Ejector Intensiv Swingle		1 pcs 27,15
L054	Spray canal cleaner needle		3 pcs 12,75
L6000	Intensiv LIPR Tray	Empty tray	1 pcs 233,00
L6000/WG69A	Intensiv LIPR Set	tray with contra-angle without light	19 pcs 1.320,00
L6000/WG69LT	Intensiv LIPR Set	tray with contra-angle with light	19 pcs 1.655,00

*prices excl VAT

Products for Stripping in Orthodontics

We thank for the clinical photographs granted:

- Dr. Giulia Fortini, Florence
- Dr. Elisabetta Carli, Florence
- Dr. Alvise Caburlotto, Venice
- Dr. Francesco Garino, Turin

Manufactured by:



Intensiv SA

Via Molino 107
6926 Montagnola
Switzerland
Tel. +41 91 986 50 50
Fax +41 91 986 50 59
info@intensiv.ch
www.intensiv.ch

Exclusively distributed in Italy by:



Leone S.p.A.

Via P. a Quaracchi 50
50019 Sesto Fiorentino
Firenze
Tel. 055 30441
Fax 055 374808
info@leone.it

Italy sales office:

Tel. 055 3044600
Fax 055 374808
italia@leone.it

www.leone.it

Products for Stripping in Orthodontics

Effective
1 february 2022

Intensiv IPR Set

Complete range
of abrasive instruments
and contra-angle



Stripping procedure

One-/double-sided interproximal reduction is a clinical procedure in use since many years. It is very useful to create space, limiting the need for extractions, in both vestibular and lingual orthodontic treatments and is necessary in aligner therapy. The Intensiv range of products dedicated to IPR (IPR=InterProximalReduction) protocols allows to treat the patient precisely, fast, and safely.

Intensiv IPR Set

- The IPR SET instruments composition is recommended and validated by Leone STEP clinical consultants. The procedure is incremental, quick and comfortable for the patient. The necessary instruments are arranged in a dedicated tray which contains the Intensiv Swingle contra-angle, a selection of Intensiv Ortho-Strips and Intensiv manual instruments, clearly arranged for intuitive use.
- The stripping protocol in orthodontics (IPR = interproximal reduction) has been developed, based on a selection of products, to perform the opening of contact points, the reduction and the finishing of proximal surfaces. The orthodontist requires efficient instruments for stripping, which are ergonomically and practically arranged for the treatment of each patient.
- The tray in polished stainless steel consists of a base and a cover connected by a hinge; its dimensions fit into the autoclaves. The closure of the tray occurs by overlapping of the cover on the base. When closed, the instruments contained in the tray are secured in place and don't move thanks to silicone supports, which facilitate the removal and organised arrangement of the instruments.
- Intensiv Ortho-Strips and Intensiv IPR-DC are inserted into two foldable internal supports, which display the article numbers of the instruments. These supports must be oriented parallel to the working surface during usage of the instruments, and are put back inside before the Tray is closed. Intensiv Ortho-Strips are arranged according to the usage sequence foreseen by the protocol marked on the tray.



DIMENSIONS
190mm length
100mm width
25m height



Content:

INTENSIV SWINGLE

ORTHO-STRIPS

MANUAL INSTRUMENTS

selected and sorted to perform precise, predictable and comfortable stripping for the patient



Intensiv ApproxOpener

Manual metal strip, serrated and with one-sided diamond coating (8µm) or double-sided diamond coating, for the opening of tight contact points. Introduction of the strip between the contact points without friction thanks to the ultra-fine grain and serrated strip border

The choice between the two instruments is based on the need for monolateral treatment, or for the treatment of both adjacent interproximal surfaces.



Intensiv manual Ortho-Strips


Diamond-coated on one side. Sterilisable, reusable.



- One-sided diamond coated strip, perforated on its working surface with two different grit sizes, 80µm (blue, coarse) and 40µm (red, medium) for efficient manual space creation.

BENEFITS

- Friction-free application due to three-hole perforation of the diamond coated strip surface
- Fast and efficient removal of excess material due to the coarse/medium grain size
- Two perforated triangular anatomical ends for a better gripping

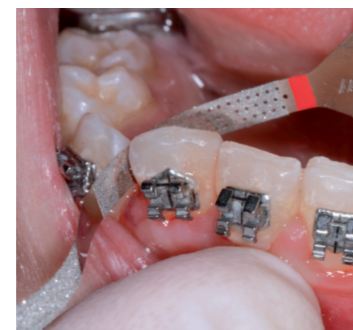

- One-sided diamond strip with two different grit sizes, 15µm (yellow, extra fine) and 8µm (orange, ultrafine), for finishing and polishing of enamel.

BENEFITS

- Efficient polishing of approximal enamel due to the extra fine and ultrafine grit size
- 2 grips for improved manual handling during treatment



Intensiv IPR-DistanceControl

Manual stainless-steel instrument to support the control of the teeth distance in the various stages of the IPR technique in orthodontics



BENEFITS

- Durable stainless steel
- Secure application due to the handle
- Variable application of the 8 instruments, individually or fixed together
- Thickness: 0,10 mm, 0,15 mm, 0,20 mm, 0,30 mm, 0,50 mm



ApproxOpener

Ortho-Strips Manual

IPR-Distance Control

IPR-DistanceControl

Support for the control of the teeth distance in the various stages of the IPR technique in orthodontics.


- REF IPRDC010/3 Thickness 0,10 mm
- REF IPRDC015/3 Thickness 0,15 mm
- REF IPRDC020/3 Thickness 0,20 mm
- REF IPRDC040/3 Thickness 0,30 mm
- REF IPRDC050/3 Thickness 0,50 mm

Intensiv Swingle Ejector
to facilitate the release of Strips from the contra-angle
REF L053

Intensiv Swingle
contra-angle with reciprocating movement
REF LWG-69 LT con luce
REF LWG-69 A senza luce

Ortho-Strips
Oscillating double-sided diamond coated strips

Ortho-Strips
One-sided diamond coated strip

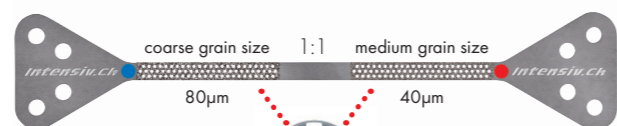

COARSE/MEDIUM
for reduction and finishing of interproximal enamel for efficient manual space creation
REF LOSM8040/6



EXTRA/ULTRAFINE
for finishing and polishing of enamel
REF LOSM1508/6


ApproxOpener
for opening of tight contact points



- 8µm Ultrafine
- one-sided diamond coating
REF LAO2018OS/3




- 8µm Ultrafine
- double-sided diamond coating
REF LAO2018DS/3

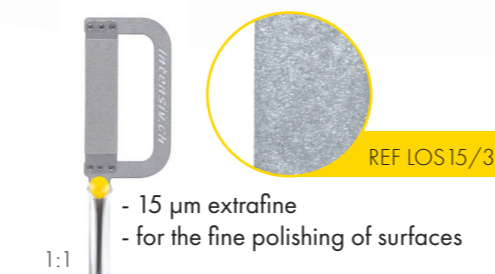

- 15 µm extrafine
- for the fine polishing of surfaces
REF LOS15/3

- 8µm Ultrafine, opener, serrated, for opening of contact points
REF LOS08OP/3
- 25 µm fine
- for step one, creating additional space
REF LOS25/3
- 60µm, medium,
- for step two, creating additional space
REF LOS60/3
- 80µm, extracoarse,
- perforated strip
- high cutting efficiency
REF LOS80C/3
- 25 µm fine
- for finishing of approximal surfaces after 60 µm and 80 µm
REF LOS25/3

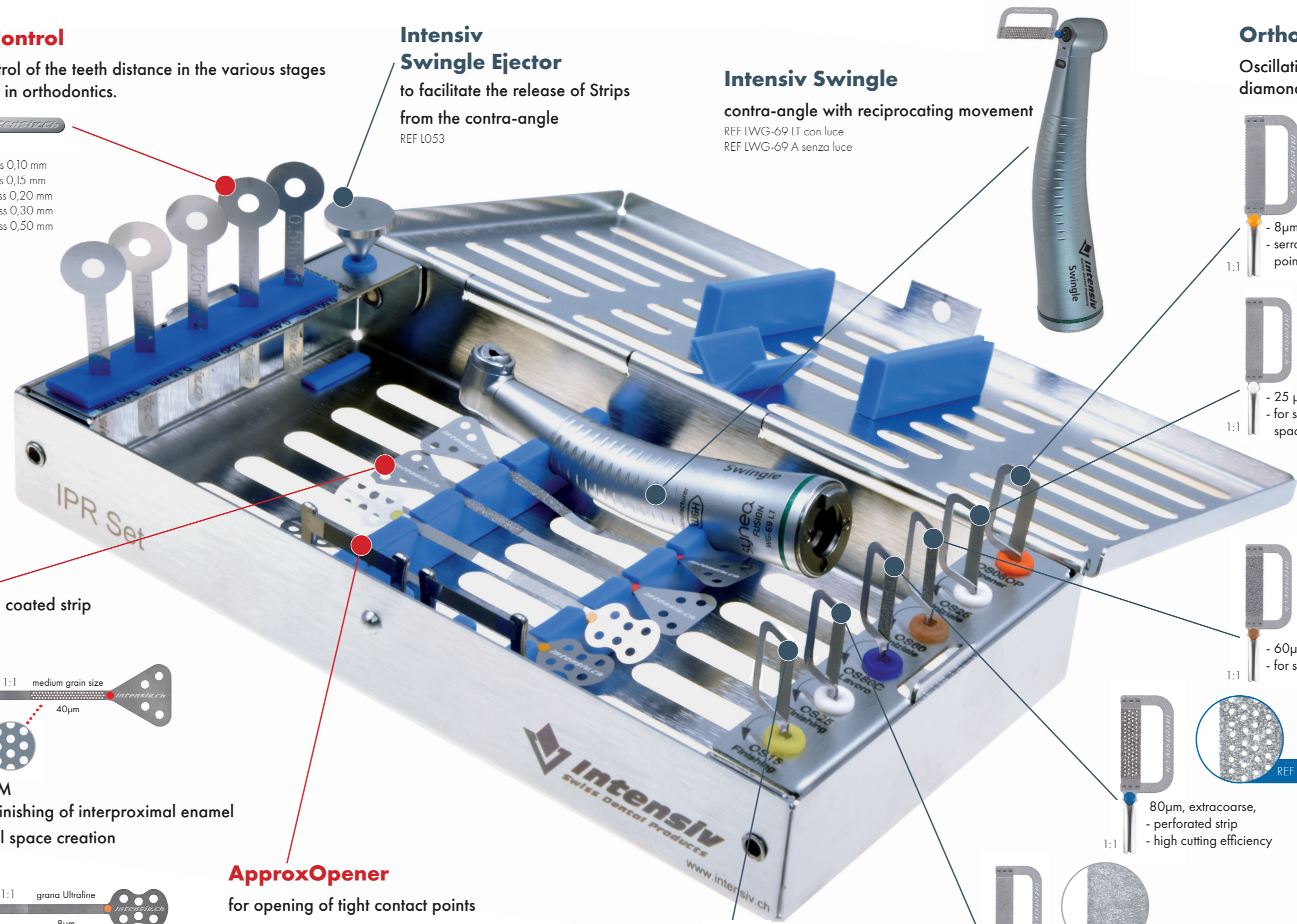

Intensiv Swingle



The reciprocating contra-angle Intensiv Swingle, manufactured by the Austrian company W&H, is, in combination with diamond coated Intensiv Ortho-Strips, the first mechanical system that efficiently meets different requirements during stripping procedures

- facilitates the rapid opening of tight interdental contact points
- enables controlled and calibrated tooth reduction
- enables precise finishing and polishing of enamel
- ensures healthy soft tissue
- safeguards patient comfort thanks to the high speed of the reciprocating movement.

PRODUCT DESCRIPTION

- Stroke 0.9 mm; 20'000 oscillating movements with 40'000 RPM
- 12 different axial fixed positions on 360°
- Contra-angle with integrated water spray and optional light
- Free axial rotating of the Ortho-Strips, to select the exact positioning of the strip

BENEFITS

- Extensive visibility on the treated area due to the slim contra-angle head and to the light
- 20.000 oscillating movements per minute for maximum performance of files and strips

Intensiv Swingle

Ejector



Intensiv Swingle Ejector

To release the Strips from the contra-angle in between the various Stripping phases

Ortho-Strips



Intensiv Ortho-Strips

Intensiv Ortho-Strips contained in the Tray are diamond coated on both sides. They are a selection of strips to be used in sequence, on the basis of clinical need for interproximal enamel reduction, therefore on the basis of the clinical space required by the method: aligners, brackets, or removable devices. Intensiv Ortho-Strips are applied for opening the interdental contact points, for reduction, and polishing of enamel.

- Refill for IPR Set composed by:
- 1 Ortho-Strips 8 µm grit Ultrafine
 - 2 Ortho-Strips 25 µm grit Fine
 - 1 Ortho-Strips 60 µm grit Medium
 - 1 Ortho-Strips 80 µm grit Coarse
 - 1 Ortho-Strips 15 µm Extrafine

