

ORTHODONTICS

2022

CLINICAL PRODUCTS



LEONEAMERICA DENTAL PRODUCTS, INC. IS LOCATED IN OXNARD, CALIFORNIA AND IS THE US ORTHODONTIC DIVISION OF LEONE ORTHODONTICS AND IMPLANTOLOGY OF FLORENCE, ITALY. LEONE WAS ESTABLISHED IN 1934 AND IS NOW ONE OF THE FEW QUALITY MANUFACTURERS PRODUCING A FULL RANGE OF ORTHODONTIC PRODUCTS. SINCE 1992, LEONE HAS BEEN A MEMBER OF THE OMA (ORTHODONTIC MANUFACTURERS ASSOCIATION) – AN ORGANIZATION FEATURING THE NINE LEADING ORTHODONTIC MANUFACTURERS WORLDWIDE. LEONEAMERICA DENTAL PRODUCTS, INC. DISTRIBUTES EXCLUSIVELY THE LEONE ORTHODONTIC RANGE IN USA, CANADA, AND PUERTO RICO.



LeoneAmerica Dental Products, Inc., Distribution site in Oxnard, CA, USA



Leone, S.p.A., Headquarters and manufacturing site in Florence, Italy



805-487-9860

Email: info@leoneamerica.com

www.leoneamerica.com

1250 Stellar Drive
Oxnard, CA 93033-2404

LeoneAmerica sole distributor of
Leone S.p.A. orthodontic products in USA,
Canada, and Puerto Rico.



CLINICAL PRODUCTS



5



21



24



37



47



50




53

D.B. BRACKETS DIAGONAL ROUND WITH FDI IDENTIFICATION

DIAGONAL ROUND ANGULATION

FOR EASY BRACKET PLACEMENT


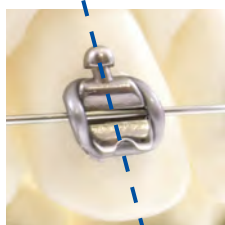


ANATOMICAL 80 MESH PAD WITH FDI IDENTIFICATION

FOR EXCELLENT BOND STRENGTH
WITH THE USE OF ANY ADHESIVE

MAXIMUM COMFORT

THANKS TO REDUCED SIZES
AND ROUNDED EDGES

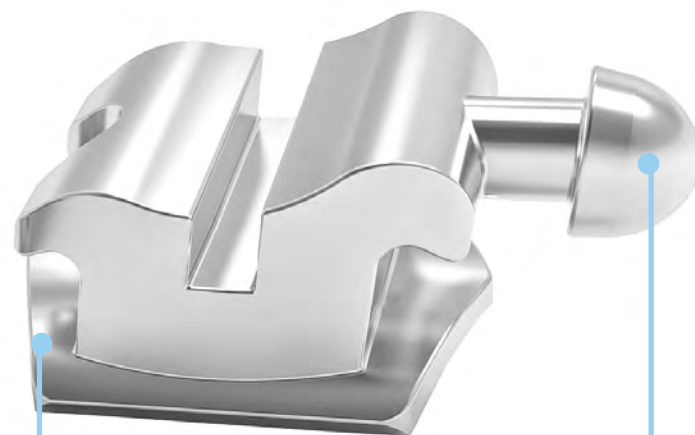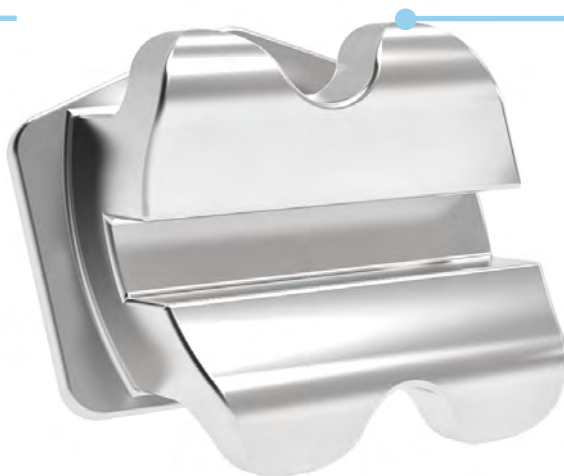

**MIM® INTEGRAL
BALL HOOK**
ON CUSPIDS
AND BICUSPIDS

































































D.B. BRACKETS DIAGONAL ROUND WITH FDI IDENTIFICATION

Manufactured using **MIM®** technology, in stainless steel welded on 80 gauge mesh sintered pads to ensure a secure bond strength to the tooth. Diagonal Round angulation, torque-in-base and calibrated sizes optimize biomechanics thus decreasing therapeutic time.

Packs of 10



CODE	CODE			pkgs.	PRICE
 .018"x.030"	 .022"x.030"		torque	ang.	
F6380-11	F6320-11	  	+12°	+5°	<u>1</u>
F6380-21	F6320-21	  	+12°	+5°	<u>1</u>
F6380-12	F6320-12	  	+8°	+9°	<u>2</u>
F6380-22	F6320-22	  	+8°	+9°	<u>2</u>
F6374-13	F6314-13	  	-2°	+13°	<u>3</u>
F6374-23	F6314-23	  	-2°	+13°	<u>3</u>
F6379-14	F6319-14	  	-7°	0°	<u>4</u>
F6379-24	F6319-24	  	-7°	0°	<u>4</u>
F6379-15	F6319-15	  	-7°	0°	<u>5</u>
F6379-25	F6319-25	  	-7°	0°	<u>5</u>
F6380-41	F6320-41	  	0°	0°	<u>1</u>
F6380-31	F6320-31	  	0°	0°	<u>1</u>
F6380-42	F6320-42	  	0°	0°	<u>2</u>
F6380-32	F6320-32	  	0°	0°	<u>2</u>
F6374-43	F6314-43	  	-11°	+7°	<u>3</u>
F6374-33	F6314-33	  	-11°	+7°	<u>3</u>
F6379-44	F6319-44	  	-17°	0°	<u>4</u>
F6379-34	F6319-34	  	-17°	0°	<u>4</u>
F6379-45	F6319-45	  	-22°	0°	<u>5</u>
F6379-35	F6319-35	  	-22°	0°	<u>5</u>



CODE	PRICE	CODE	PRICE	
20 brackets 1 case		200 brackets 10 case		inch
F6380-91	\$65.50	F6381-91	\$491.40	.018"
F6320-91	\$65.50	F6321-91	\$491.40	.022"

SPECIAL
BUY 2 OR MORE
MIX AND MATCH
FOR **\$30.00**
EACH

SPECIAL:
BUY 2
CASES FOR
\$60.00
PER CASE


D.B. METAL BRACKETS WITH FDI IDENTIFICATION

ANATOMICAL 80 MESH PAD

FOR EXCELLENT BOND STRENGTH WITH THE USE OF ANY ADHESIVE, BICUSPIDS FEATURE A BONDING PAD EXTENDED TO OCCLUSAL TO MINIMIZE INTERFERENCES

BALL HOOK

MIM® INJECTION MOLDED AND INTEGRATED ON CUSPIDS AND BICUSPIDS

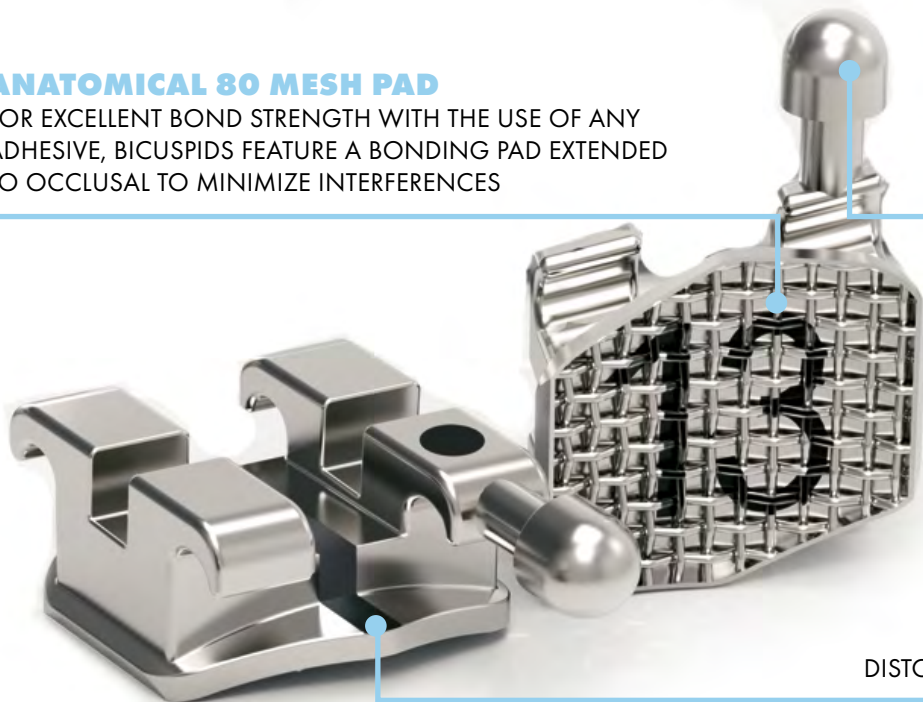


CENTRAL SCRIBE LINE

LASER MARKED BETWEEN TIE WINGS AND IDENTIFICATION
DISTO-GINGIVAL DOT FOR EASE OF PLACEMENT



FDI SYSTEM MARKING ON THE MESH PAD

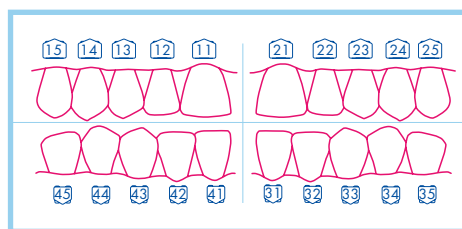


















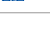










































An FDI (Fédération Dentaire Internationale) identification number is laser marked on the mesh pad of every bracket, clearly indicating the tooth to which a specific bracket will be attached. This method identifies each tooth with a two-digit Arabic number: the first digit indicates the quadrant, the second identifies the tooth according to Palmer's notation. For example: the bracket for the upper right central incisor has number 11 laser-marked on the mesh pad. The number, permanently etched and free of toxic substances, is immediately recognizable and easily readable without the aid of magnifying lenses.

The FDI system laser marking also improves bond strength.

MINI DIAGONAL® ROTH SYSTEM



CODE	CODE			torque	ang.		pkg.	PRICE
 .018"x.030"	 .022"x.030"							
F7281-02	F7221-02	  3.3		+12°	+5°	<u>1</u>	10	\$37.20
F7281-03	F7221-03	  3.3		+12°	+5°	<u>1</u>	10	\$37.20
F7282-02	F7222-02	  2.7		+8°	+9°	<u>2</u>	10	\$37.20
F7282-03	F7222-03	  2.7		+8°	+9°	<u>2</u>	10	\$37.20
F7283-02	F7223-02	  2.7		-2°	+13°	<u>3</u>	10	\$37.20
F7283-03	F7223-03	  2.7		-2°	+13°	<u>3</u>	10	\$37.20
F7273-02	F7213-02	  2.7		-2°	+13°	<u>3</u>	10	\$37.20
F7273-03	F7213-03	  2.7		-2°	+13°	<u>3</u>	10	\$37.20
F7184-02	F7124-02	  2.7		-7°	0°	<u>4</u>	10	\$37.20
F7184-03	F7124-03	  2.7		-7°	0°	<u>4</u>	10	\$37.20
F7274-02	F7214-02	  2.7		-7°	0°	<u>4</u>	10	\$37.20
F7274-03	F7214-03	  2.7		-7°	0°	<u>4</u>	10	\$37.20
F7185-02	F7125-02	  2.7		-7°	0°	<u>5</u>	10	\$37.20
F7185-03	F7125-03	  2.7		-7°	0°	<u>5</u>	10	\$37.20
F7179-02	F7119-02	  2.7		-7°	0°	<u>5</u>	10	\$37.20
F7179-03	F7119-03	  2.7		-7°	0°	<u>5</u>	10	\$37.20
F7081-06	F7021-06	  2.1		0°	0°	<u>1</u>	10	\$37.20
F7081-07	F7021-07	  2.1		0°	0°	<u>1</u>	10	\$37.20
F7082-06	F7022-06	  2.1		0°	0°	<u>2</u>	10	\$37.20
F7082-07	F7022-07	  2.1		0°	0°	<u>2</u>	10	\$37.20
F7283-06	F7223-06	  2.7		-11°	+7°	<u>3</u>	10	\$37.20
F7283-07	F7223-07	  2.7		-11°	+7°	<u>3</u>	10	\$37.20
F7273-06	F7213-06	  2.7		-11°	+7°	<u>3</u>	10	\$37.20
F7273-07	F7213-07	  2.7		-11°	+7°	<u>3</u>	10	\$37.20
F7184-06	F7124-06	  2.7		-17°	0°	<u>4</u>	10	\$37.20
F7184-07	F7124-07	  2.7		-17°	0°	<u>4</u>	10	\$37.20
F7274-06	F7214-06	  2.7		-17°	0°	<u>4</u>	10	\$37.20
F7274-07	F7214-07	  2.7		-17°	0°	<u>4</u>	10	\$37.20
F7185-06	F7125-06	  2.7		-22°	0°	<u>5</u>	10	\$37.20
F7185-07	F7125-07	  2.7		-22°	0°	<u>5</u>	10	\$37.20
F7175-06	F7115-06	  2.7		-22°	0°	<u>5</u>	10	\$37.20
F7175-07	F7115-07	  2.7		-22°	0°	<u>5</u>	10	\$37.20




Brackets
not available in kits

SPECIAL
BUY 2 OR MORE
MIX AND MATCH
FOR **\$34.00**
EACH

MINI DIAGONAL® ROTH SYSTEM













CODE	PRICE	CODE	PRICE	CODE	PRICE	CODE	PRICE	
20 brackets 1 case		200 brackets 10 case		500 brackets 25 cases		1000 brackets 50 cases		
F7270-91	\$78.10	F7271-91	\$741.50	F7273-91	\$1,446.90			.018"
F7210-91	\$78.10	F7211-91	\$741.50	F7213-91	\$1,446.90	F7215-91	\$2,730.00	.022"

**BUY 2
OR MORE
CASES
GET 20%
OFF**

MINI DIAGONAL® MBT* SYSTEM



CODE		CODE				pkg.		PRICE	
 .018"x.030"	 .022"x.030"			torque	ang.				
F7380-11	F7320-11				+17°	+4°	<u>1</u>	10	\$37.20
F7380-21	F7320-21				+17°	+4°	<u>1</u>	10	\$37.20
F7380-12	F7320-12				+10°	+8°	<u>2</u>	10	\$37.20
F7380-22	F7320-22				+10°	+8°	<u>2</u>	10	\$37.20
F7380-13	F7320-13				-7°	+8°	<u>3</u>	10	\$37.20
F7380-23	F7320-23				-7°	+8°	<u>3</u>	10	\$37.20
F7370-13	F7310-13				-7°	+8°	<u>3</u>	10	\$37.20
F7370-23	F7310-23				-7°	+8°	<u>3</u>	10	\$37.20
F7380-14	F7320-14				-7°	0°	<u>4</u>	10	\$37.20
F7380-24	F7320-24				-7°	0°	<u>4</u>	10	\$37.20
F7370-14	F7310-14				-7°	0°	<u>4</u>	10	\$37.20
F7370-24	F7310-24				-7°	0°	<u>4</u>	10	\$37.20
F7380-15	F7320-15				-7°	0°	<u>5</u>	10	\$37.20
F7380-25	F7320-25				-7°	0°	<u>5</u>	10	\$37.20
F7370-15	F7310-15				-7°	0°	<u>5</u>	10	\$37.20
F7370-25	F7310-25				-7°	0°	<u>5</u>	10	\$37.20
F7380-41	F7320-41				-6°	0°	<u>1</u>	10	\$37.20
F7380-31	F7320-31				-6°	0°	<u>1</u>	10	\$37.20
F7380-42	F7320-42				-6°	0°	<u>2</u>	10	\$37.20
F7380-32	F7320-32				-6°	0°	<u>2</u>	10	\$37.20
F7370-43	F7320-43				-6°	+3°	<u>3</u>	10	\$37.20
F7370-33	F7320-33				-6°	+3°	<u>3</u>	10	\$37.20
F7370-43	F7310-43				-6°	+3°	<u>3</u>	10	\$37.20
F7370-33	F7310-33				-6°	+3°	<u>3</u>	10	\$37.20
F7380-44	F7320-44				-12°	+2°	<u>4</u>	10	\$37.20
F7380-34	F7320-34				-12°	+2°	<u>4</u>	10	\$37.20
F7370-44	F7310-44				-12°	+2°	<u>4</u>	10	\$37.20
F7370-34	F7310-34				-12°	+2°	<u>4</u>	10	\$37.20
F7380-45	F7320-45				-17°	+2°	<u>5</u>	10	\$37.20
F7380-35	F7320-35				-17°	+2°	<u>5</u>	10	\$37.20
F7370-45	F7310-45				-17°	+2°	<u>5</u>	10	\$37.20
F7370-35	F7310-35				-17°	+2°	<u>5</u>	10	\$37.20




Brackets
not available in kits

SPECIAL
BUY 2 OR MORE
MIX AND MATCH
FOR **\$32.00**
EACH

MINI DIAGONAL® MBT* SYSTEM



CODE	PRICE	CODE	PRICE	CODE	PRICE	CODE	PRICE	
20 brackets 1 case		200 brackets 10 case		500 brackets 25 cases		1000 brackets 50 cases		
F7370-91	\$78.10	F7371-91	\$741.50	F7373-91	\$1,446.90	F7375-91	\$2,730.00	.018"
F7310-91	\$78.10	F7311-91	\$741.50	F7313-91	\$1,446.90	F7315-91	\$2,730.00	.022"

**BUY 2
OR MORE
CASES
GET 20%
OFF**

D.B. BRACKETS EXTREMO NO-NICKEL

EXTREMO

Single piece bracket **MIM®**

- **PRECISION**
- **CONSTANT QUALITY**
- **NO SOLDERING POINTS**
- **COMPUTER DESIGNED**
- **ANATOMICAL BOND BASES**
- EXCELLENT BOND STRENGTH**
- **BIOCOMPATIBILITY**
- **PATIENT COMFORT**





FDI SYSTEM MARKING
ON THE MESH PAD
(FOR ROTH PRESCRIPTION)

HIGH RETENTION BASE
WITH MICRO CAVITIES
TO ENSURE STRONG RETENTION TO
THE TOOTH



MIM® TECHNOLOGY
PERFECTLY ROUNDED EDGES
FOR MAXIMUM COMFORT

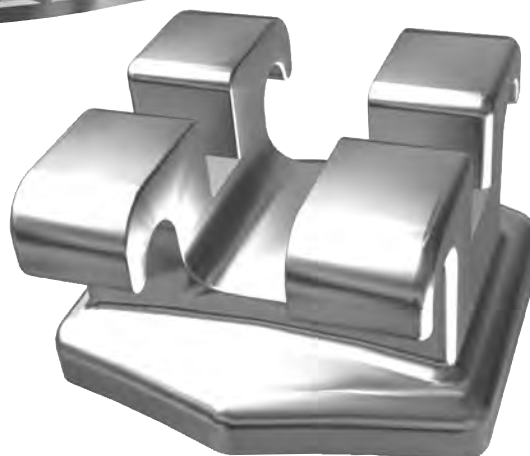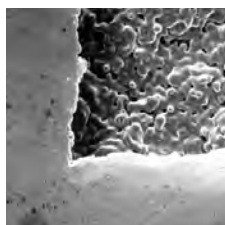










































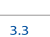





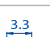











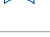


Extremo nickel-free* brackets are entirely biocompatible. They are CAD designed to get anatomical fits and **MIM®** moulded in single piece, without any soldering point. The anatomical base features highly retentive micro cavities for a perfect bond strength and rounded edges to provide patient compliance. The base on Roth prescription brackets is laser marked with FDI number for immediate identification.

*Nickel trace elements of $\leq 0,1\%$ are due to the process of manufacturing

EXTREMO NO-NICKEL ROTH SYSTEM WITH FDI IDENTIFICATION



CODE	CODE		torque	ang.		pkg.	PRICE
 .018"x.030"	 .022"x.030"	.					
F9280-11	F9220-11	   3.6	+12°	+5°	<u>1</u>	10	\$31.50
F9280-21	F9220-21	   3.6	+12°	+5°	<u>1</u>	10	\$31.50
F9280-12	F9220-12	   3	+8°	+9°	<u>2</u>	10	\$31.50
F9280-22	F9220-22	   3	+8°	+9°	<u>2</u>	10	\$31.50
F9270-13	F9210-13	   3.3	-2°	+13°	<u>3</u>	10	\$31.50
F9270-23	F9210-23	   3.3	-2°	+13°	<u>3</u>	10	\$31.50
F9270-14	F9210-14	   3.3	-7°	0°	<u>4</u>	10	\$31.50
F9270-24	F9210-24	   3.3	-7°	0°	<u>4</u>	10	\$31.50
F9270-15	F9210-15	   3.3	-7°	0°	<u>5</u>	10	\$31.50
F9270-25	F9210-25	   3.3	-7°	0°	<u>5</u>	10	\$31.50
F9280-41	F9220-41	   2.6	0°	0°	<u>1</u>	10	\$31.50
F9280-31	F9220-31	   2.6	0°	0°	<u>1</u>	10	\$31.50
F9280-42	F9220-42	   2.6	0°	0°	<u>2</u>	10	\$31.50
F9280-32	F9220-32	   2.6	0°	0°	<u>2</u>	10	\$31.50
F9270-43	F9210-43	   3.3	-11°	+7°	<u>3</u>	10	\$31.50
F9270-33	F9210-33	   3.3	-11°	+7°	<u>3</u>	10	\$31.50
F9270-44	F9210-44	   3.3	-17°	0°	<u>4</u>	10	\$31.50
F9270-34	F9210-34	   3.3	-17°	0°	<u>4</u>	10	\$31.50
F9270-45	F9210-45	   3.3	-22°	0°	<u>5</u>	10	\$31.50
F9270-35	F9210-35	   3.3	-22°	0°	<u>5</u>	10	\$31.50

CODE	PRICE	CODE	PRICE	CODE	PRICE	CODE	PRICE	
20 brackets 1 case		200 brackets 10 case		500 brackets 25 cases		1000 brackets 50 cases		
F9270-91	\$62.70	F9271-91	\$597.50	F9273-91	\$1,448.20	F9275-91	\$2,835.00	.018"
F9210-91	\$62.70	F9211-91	\$597.50	F9213-91	\$1,448.20	F9215-91	\$2,835.00	.022"

SPECIAL
BUY 2 OR MORE
MIX AND MATCH
FOR **\$28.00**
EACH

SPECIAL
BUY 2 CASES
FOR
\$49.00
EACH

D.B. BUTTONS EXTREMO NO-NICKEL

Single piece buttons manufactured with **MIM®** technology in a biocompatible **nickel free*** alloy. The Extremo no-nickel base, anatomical and with micro highly retentive cavities, perfectly fits the tooth and ensures an excellent bond with the use of any kind of adhesives. Packs of 10



F9060-00
Low flat
PRICE \$26.00



F9070-00
Deep curved
PRICE \$26.00


SPECIAL
BUY 2
FOR
\$18.00
PER PACK

D.B. TUBES EXTREMO NO-NICKEL

Single piece tubes manufactured with **MIM®** technology in a biocompatible nickel free* alloy. The Extremo No-Nickel base, anatomical with micro highly retentive cavities, perfectly fits the tooth and ensures an excellent bond with the use of any kind of adhesives. Packs of 10 - Packs of 100 on request


FUNNELED ENTRY WAY
FOR EASE OF WIRE INSERTION



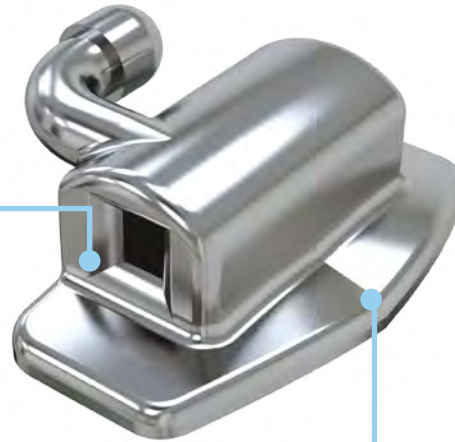





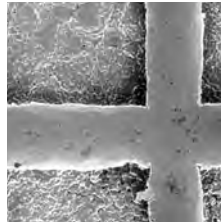














ANATOMICAL MESH BASE
WITH MICRO CELL
INTERNAL SURFACE



LOW PROFILE DESIGN
THAT HELPS REDUCE
OCCLUSAL INTERFERENCE

OCCLUSAL INDENT
FOR A PERFECT ANATOMICAL
FIT TO MOLARS




CODE								pkg.	PRICE		
						tor- que	rotat.	slot	direct bond- ing tubes		
	F8581-12	$\frac{7}{7} \frac{6}{6} \frac{7}{7}$				0°	0°	.018"	R	10	\$21.50
	F8581-13					0°	0°	.018"	L	10	\$21.50
	F8521-12					0°	0°	.022"	R	10	\$21.50
	F8521-13					0°	0°	.022"	L	10	\$21.50
	F8681-12	$\frac{7}{7} \frac{6}{6} \frac{7}{7}$				-10°	+8°	.018"	R	10	\$21.50
	F8681-13					-10°	+8°	.018"	L	10	\$21.50
	F8621-22					-10°	+8°	.022"	R	10	\$21.50
	F8621-23					-10°	+8°	.022"	L	10	\$21.50
	F8621-32	$\frac{7}{7} \frac{6}{6} \frac{7}{7}$				-14°	+8°	.022"	R	10	\$21.50
	F8621-33					-14°	+8°	.022"	L	10	\$21.50
	F8621-26	$\frac{7}{7} \frac{6}{6} \frac{7}{7}$				-20°	+8°	.022"	R		\$21.50
	F8621-27					-20°	+8°	.022"	L	10	\$21.50
	F8681-16	$\frac{7}{7} \frac{6}{6} \frac{7}{7}$				-25°	+8°	.018"	R	10	\$21.50
	F8681-17					-25°	+8°	.018"	L	10	\$21.50
	F8621-36					-25°	+8°	.022"	R	10	\$21.50
	F8621-37					-25°	+8°	.022"	L	10	\$21.50
	F8570-02	$\frac{6}{6} \frac{6}{6}$				0°	0°	.018"	R	10	\$71.20
	F8570-03					0°	0°	.018"	L	10	\$71.20
	F8510-02	$\frac{6}{6} \frac{6}{6}$				0°	0°	.022"	R	10	\$71.20
	F8510-03					0°	0°	.022"	L	10	\$71.20

**SPECIAL
BUY 2
FOR
\$44.00**

SLIDE™ LOW FRICTION LIGATURES





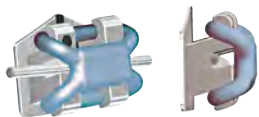




SLIDE is a product manufactured with a special polyurethane mix for medical use. It is applied similarly to conventional elastic ligatures. Once placed on the bracket, it creates a passive ligature on the slot that makes the archwire free to **SLIDE** and to produce its effects on the dentoalveolar structures. This ligature can be used to create low friction with any bracket, but is especially recommended with Logic, Mini and Midi Diagonali, Aqua and Idea-L bracket lines. The particular conformation of the **SLIDE** ligature considerably improves patient's comfort during the first phases of treatment.













Packs of 10 modules, 6 ligatures each.



LEARN
MORE



CODE	CODE	CODE	colour		PRICE
extra short	short	medium			
K6254-10A	K6251-10A	K6252-10A		blue	\$80.60
K6254-10G	K6251-10G	K6252-10G		yellow	\$80.60
K6254-10V	K6251-10V	K6252-10V		green	\$80.60
K6254-10	K6251-10	K6252-10		white	\$80.60
K6254-10R	K6251-10R	K6252-10R		red	\$80.60
K6254-10S	K6251-10S	K6252-10S		pink	\$80.60
K6254-10N	K6251-10N	K6252-10N		black	\$80.60
K6224-10	K6221-10	K6222-10		ice	\$80.60
K6264-10	K6261-10	K6262-10		silver	\$80.60

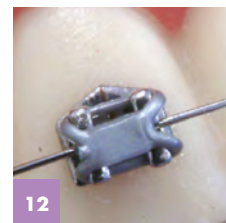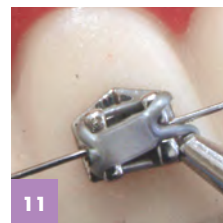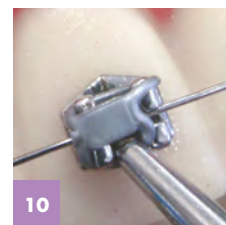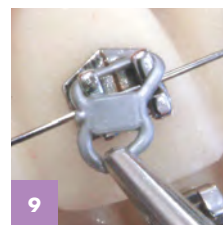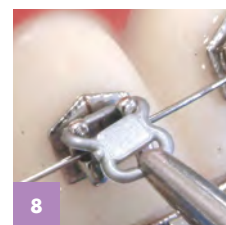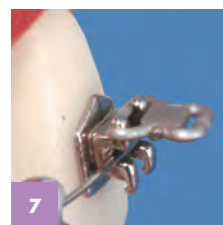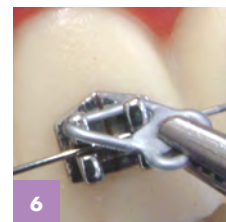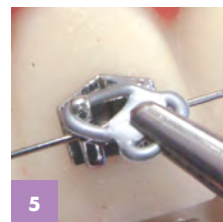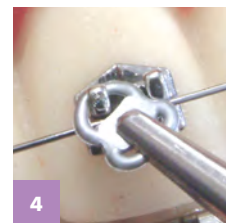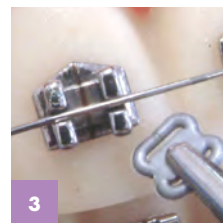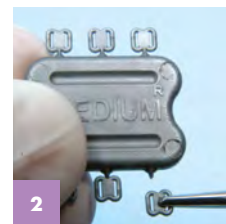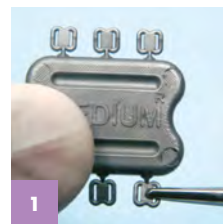

INTERACTIVE D.B. METAL BRACKETS SELF-LIGATING




Inter.\activeSL

BALL HOOKS:
PERFECTLY SPHERICAL
AND LOW PROFILE
FOR MAXIMUM COMFORT.




TIE WINGS:
FOR THE APPLICATION
OF ANY KIND OF LIGATURE
IN THE EVENTUAL NEED FOR
A BIOMECHANICAL TOTAL
CONTROL.

SPRING CLIP:
MADE OF A HIGHLY ELASTIC ALLOY,
IT GUARANTEES EASY AND SAFE
OPEN/CLOSE MOVEMENTS.
THE SLIDING PORTION OF THE CLIP
IS CONSTRAINED BETWEEN THE BASE
AND THE BODY OF THE BRACKET BY
A LASER WELDING THAT PREVENTS
THE POSSIBILITY OF DETACHMENT

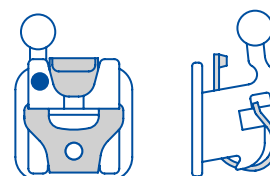

**BEVELED SLOT
ON MESIO-DISTAL EDGE:**
ALLOWS SLIDING OF ARCHWIRES
FOR HEAVY MISALIGNMENTS,
WHILE AVOIDING NOTCHING
AND BINDING.

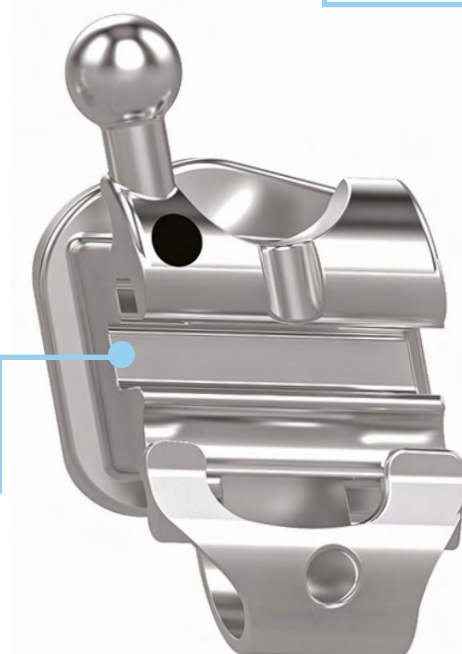

MIM® TECHNOLOGY:
THE IDEAL TECHNOLOGY FOR
THE COMPLEXITY OF THE DESIGN,
MAXIMUM ACCURACY OF THE
SLOT AND UNDER THE TIE WINGS.

MAXIMUM COMFORT AND VERSATILITY

InterActive SL self-ligating brackets, due to their design and reduced dimensions, offer a very comfortable treatment for the patient without loss of control. The interactive clip ensures an efficient clinical management with progressive biomechanical control to take advantage of the low friction early in treatment stages, while achieving perfect finishing with the final archwires.

MATERIAL & DESIGN

The bracket body is made of biomedical steel and is laser welded on the 80 mesh gauge base with FDI identification. The colored dot with disto-gingival identification indicates each bracket quadrant, with the same color code of all Leone bracket range.

INTERACTIVE CONTROL

The specific shape of clip and slot of InterActive brackets allows the user to modulate the most appropriate level of friction force between bracket and wire, depending on the needs of the various stages of treatment.



PASSIVE PHASE

First stage round arches are not bound by the clip: the low friction will facilitate the process of alignment and leveling



INTERACTIVE PHASE

Rectangular arch wires, used for space closure, rotation, and torque control, work to elastically deform the clip for the biomechanical control necessary in these stages of treatment



ACTIVE PHASE

Rectangular arches for finishing and detailing completely fill the slot while going into active contact with clip: this allows the exploitation of elastic properties in order to obtain minimal movements for finishing of treatment.


EASY OPEN/CLOSE

The clip has a central hole and does not require any special tool for opening and closing.



OPENING

Insert the tip of a probe or utility tool into the hole in the clip and move towards the occlusal plane

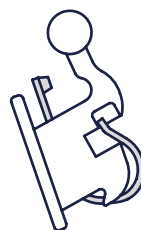

CLOSING






























































Slide the clip with a slight pressure towards the gums using a tool tip or even just a finger

INTERACTIVE D.B. METAL BRACKETS SELF-LIGATING



INTERACTIVE SELF LIGATING D.B. BRACKETS ROTH SYSTEM

**SPECIAL
BUY 4
GET 1
FREE**

CODE				torque	ang.		pkg.	PRICE
 .022"x .030"								
F1100-11				+12°	+5°	1	5	\$36.60
F1100-21				+12°	+5°	1	5	\$36.60
F1100-12				+8°	+9°	2	5	\$36.60
F1100-22				+8°	+9°	2	5	\$36.60
F1100-13				-2°	+13°	3	5	\$36.60
F1100-23				-2°	+13°	3	5	\$36.60
F1100-14				-7°	0°	4	5	\$36.60
F1100-24				-7°	0°	4	5	\$36.60
F1100-15				-7°	0°	5	5	\$36.60
F1100-25				-7°	0°	5	5	\$36.60
F1100-41				0°	0°	1	5	\$36.60
F1100-31				0°	0°	1	5	\$36.60
F1100-42				0°	0°	2	5	\$36.60
F1100-32				0°	0°	2	5	\$36.60
F1100-43				-11°	+7°	3	5	\$36.60
F1100-33				-11°	+7°	3	5	\$36.60
F1100-44				-17°	0°	4	5	\$36.60
F1100-34				-17°	0°	4	5	\$36.60
F1100-45				-22°	0°	5	5	\$36.60
F1100-35				-22°	0°	5	5	\$36.60




































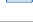

























INTERACTIVE SELF LIGATING D.B. BRACKETS KIT ROTH SYSTEM



CODE		PRICE	CODE		PRICE
1 case 20 Brackets			10 kits of 1 case		
F1100-91	.022"x .030"	\$146.00	F1101-91	.022"x .030"	\$1,167.60

INTERACTIVE SELF LIGATING D.B. BRACKETS MBT* SYSTEM

**SPECIAL
BUY 4
GET 1
FREE**


CODE				torque	ang.		pkg.	PRICE
 .022"x .030"								
F1102-11				+17°	+4°	1	5	\$36.60
F1102-21				+17°	+4°	1	5	\$36.60
F1102-12				+10°	+8°	2	5	\$36.60
F1102-22				+10°	+8°	2	5	\$36.60
F1102-13				-7°	+8°	3	5	\$36.60
F1102-23				-7°	+8°	3	5	\$36.60
F1100-14				-7°	0°	4	5	\$36.60
F1100-24				-7°	0°	4	5	\$36.60
F1100-15				-7°	0°	5	5	\$36.60
F1100-25				-7°	0°	5	5	\$36.60
F1102-41				-6°	0°	1	5	\$36.60
F1102-31				-6°	0°	1	5	\$36.60
F1102-42				-6°	0°	2	5	\$36.60
F1102-32				-6°	0°	2	5	\$36.60
F1102-43				-6°	+3°	3	5	\$36.60
F1102-33				-6°	+3°	3	5	\$36.60
F1102-44				-12°	+2°	4	5	\$36.60
F1102-34				-12°	+2°	4	5	\$36.60
F1102-45				-17°	+2°	5	5	\$36.60
F1102-35				-17°	+2°	5	5	\$36.60

INTERACTIVE SELF LIGATING D.B. BRACKETS KIT MBT* SYSTEM



CODE		PRICE	CODE		PRICE
1 case 20 Brackets			10 kits of 1 case		
F1100-92	.022"x .030"	\$146.00	F1101-92	.022"x .030"	\$1,167.60

D.B. AQUA CERAMIC BRACKETS

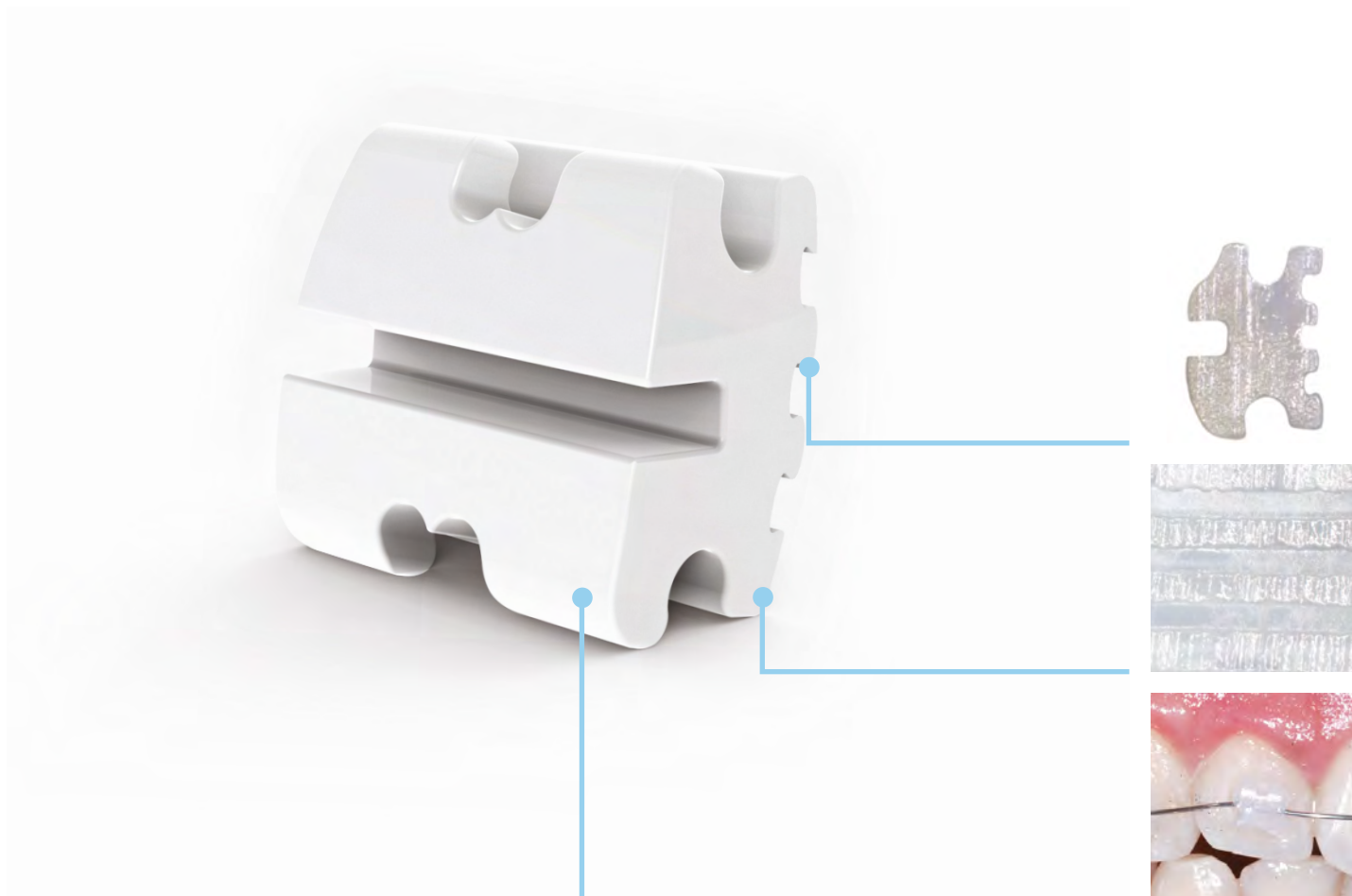

PURE AESTHETICS































































- ENTIRELY MADE USING A SINTERED ANTIALLERGENIC POLYCRYSTALLINE ALUMINA
- MAXIMUM TRANSLUCENCE AND METAL SLOT FREE FOR NATURAL AESTHETICS AGAINST ANY SHADE OF ENAMEL
- LOW FRICTION BIOMECHANICS, USED IN CONJUNCTION WITH **SLIDE™** LIGATURES
- OPTIMAL MECHANIC RETENTION WITH THE USE OF ANY DESIRED ADHESIVE
- CUSPID AND BICUSPID BRACKETS WITH BIDIRECTIONAL HOOK

D.B. AQUA CERAMIC BRACKETS


AQUA ceramic brackets are manufactured with a sintering process to achieve the purest alumina and a polycrystalline structure with an average uniform grain. The result is a bracket combining the best translucence with a superior mechanical resistance thus assuring a naturally aesthetic look to match every tooth shade. AQUA brackets are free of any metal slot and totally anti-allergenic. No specific primers are required since the base has got a new design providing a special mechanical retention, thus allowing the use of any orthodontic adhesive available. Polished surface, smooth edges and low profile design make the AQUA ceramic brackets superior in comfort.

AQUA ROTH SYSTEM



CODE	CODE		torque	ang.		pkg.	PRICE
							
.018"x.030"	.022"x.030"						
F5781-02	F5721-02	  	+12°	+5°	<u>1</u>	1	\$16.30
F5781-03	F5721-03	  	+12°	+5°	<u>1</u>	1	\$16.30
F5782-02	F5722-02	  	+8°	+9°	<u>2</u>	1	\$16.30
F5782-03	F5722-03	  	+8°	+9°	<u>2</u>	1	\$16.30
F5773-02	F5713-02	  	-2°	+13°	<u>3</u>	1	\$16.30
F5773-03	F5713-03	  	-2°	+13°	<u>3</u>	1	\$16.30
F5779-02	F5719-02	  	-7°	0°	<u>4</u>	1	\$16.30
F5779-03	F5719-03	  	-7°	0°	<u>4</u>	1	\$16.30
F5779-02	F5719-02	  	-7°	0°	<u>5</u>	1	\$16.30
F5779-03	F5719-03	  	-7°	0°	<u>5</u>	1	\$16.30
F5780-05	F5720-05	  	0°	0°	<u>1</u>	1	\$16.30
F5780-05	F5720-05	  	0°	0°	<u>1</u>	1	\$16.30
F5780-05	F5720-05	  	0°	0°	<u>2</u>	1	\$16.30
F5780-05	F5720-05	  	0°	0°	<u>2</u>	1	\$16.30
F5773-06	F5713-06	  	-11°	+7°	<u>3</u>	1	\$16.30
F5773-07	F5713-07	  	-11°	+7°	<u>3</u>	1	\$16.30
F5774-06	F5714-06	  	-17°	0°	<u>4</u>	1	\$16.30
F5774-07	F5714-07	  	-17°	0°	<u>4</u>	1	\$16.30
F5775-06	F5715-06	  	-22°	0°	<u>5</u>	1	\$16.30
F5775-07	F5715-07	  	-22°	0°	<u>5</u>	1	\$16.30



CODE	PRICE	CODE	PRICE	CODE	PRICE	CODE	PRICE	
6 brackets 3 2 1 1 2 3		10 brackets 5 4 3 2 1 1 2 3 4 5		12 brackets 3 2 1 1 2 3 3 2 1 1 2 3		20 brackets 5 4 3 2 1 1 2 3 4 5 5 4 3 2 1 1 2 3 4 5		
F5706-80	\$93.20	F5710-80	\$154.90	F5712-80	\$185.90	F5720-80	\$309.80	.018"
F5706-20	\$93.20	F5710-20	\$154.90	F5712-20	\$185.90	F5720-20	\$309.80	.022"

**BUY TWO
OR MORE
CASES
GET 20%
OFF**

IDEA-L LINGUAL BRACKETS

IDEA-L lingual brackets feature small mesio-distal sizes, 1.4 mm depth and a special bracket design which make the application of Slide ligatures easy. In most cases, the use of a round wire **MEMORIA®** .012 is sufficient to get the solution of overcrowding and the aesthetic improvement of smile. In cases of dental rotations and/or major misplacements, wires of higher diameter, like **MEMORIA®** or Beta **MEMORIA®**, and conventional ligatures may also be used to gain biomechanical control. To get perfect alignment and levelling of the anterior teeth, an accurate position of the bracket is necessary; it is therefore advisable to use an indirect bonding technique with arch transfer trays made in the laboratory or with the special "Jigs" conceived for this technique.

ideal
O LIGHT LINGUAL SYSTEM

PROFILE OF 1,4 MM DEPTH

FOR PATIENT ACCEPTANCE AND
MINIMAL SPEECH ISSUES




CUSTOM DESIGN

MAKES THE APPLICATION OF THE
SLIDE™ LIGATURE EASY AND KEEPS
IT IN PLACE DURING TREATMENT



SLIDE™ LIGATURE
EMBODIES THE PROFILE
OF THE BRACKET,
MAKING IT SMOOTH
AND COMFORTABLE

SPECIAL PLASTIC JIG (FOR RECOMMENDED INDIRECT BONDING)

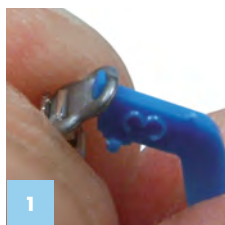

1
Assemble jig on the
bracket with the
SLIDE™ ligature
hooked under
the gingival wing




2
Measure of the height of the selected jig on working cast and select correct size for the best bracket position



3
Position in the mouth and polymerization of the adhesive





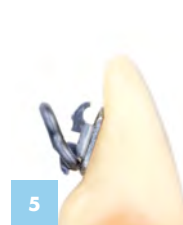

























4
Removal of the jig with a rotary movement



5
Bracket with the
SLIDE™ ligature
under the gingival wing



6
Bracket with wire and the
SLIDE™ ligature in place


CODE						pkg.	PRICE
			.018"x.030"	torque	ang.		
F4980-04				0°	0°	$\frac{1}{1}$	5 \$93.50
F4980-04				0°	0°	$\frac{2}{2}$	5 \$93.50
F4980-04				0°	0°	$\frac{3}{3}$	5 \$93.50
F4984-04				0°	0°	$\frac{4}{4}$	5 \$93.50
F4980-04				0°	0°	$\frac{1}{1}$	5 \$93.50
F4980-04				0°	0°	$\frac{2}{2}$	5 \$93.50
F4980-04				0°	0°	$\frac{3}{3}$	5 \$93.50
F4984-04				0°	0°	$\frac{4}{4}$	5 \$93.50

SPECIAL:
BUY 2 FOR
\$89.00 EACH

MEMORIA® AND BETA MEMORIA® LINGUAL ARCHWIRES

Sections .012" and .014" are made of a special Nickel Titanium alloy allowing light adaptations necessary in case of treatment extension to first bicuspids, thus keeping high elasticity. Section .016" Beta **MEMORIA®** is ideal to optimize the alignment and levelling values. The selector C4950-00 is available to select the correct archwire shape.

Packs of 2

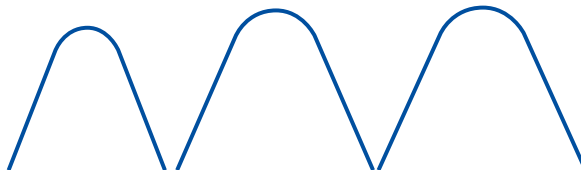

MEMORIA® LINGUAL ARCHWIRES	CODE	CODE	CODE	Ø	PRICE
	small	medium	large	inch	
	C4950-12	C4951-12	C4952-12	.012	\$12.60
	C4950-14	C4951-14	C4952-14	.014	\$12.60
BETA MEMORIA® LINGUAL ARCHWIRES	CODE	CODE	CODE	Ø	PRICE
	C4970-16	C4971-16	C4972-16	.016	\$12.60

F4980-91 IDEA-L LINGUAL BRACKET KIT

Composition of the kit: 8 brackets (one arch) and 12 **SLIDE™** ligatures extra-small, silver colour

PRICE \$89.30



F4900-60 KIT OF ASSORTED JIGS

Six sizes for the selection of the proper position according to the anatomy of the lingual surface of the tooth. Assorted pack of 120 pcs. **PRICE \$31.50**




JIGS - REFILL	CODE	size	pkg.	PRICE
	F4900-01	size 1	20	\$8.40
	F4900-02	size 2	20	\$8.40
	F4900-03	size 3	20	\$8.40
	F4900-04	size 4	20	\$8.40
	F4900-05	size 5	20	\$8.40
	F4900-06	size 6	20	\$8.40

C4950-00 LINGUAL ARCHWIRE SELECTOR

Lay the selector on the model and choose the correct arch shape.

PRICE \$8.40

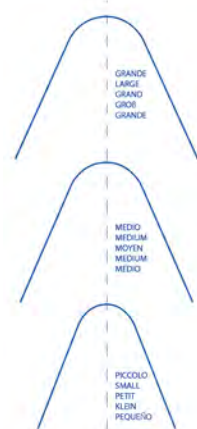

Courtesy of Dr. F. Giuntoli

D.B. MIM® SINGLE TUBES LOW PROFILE WITH FDI IDENTIFICATION

These D.B. single tubes, made using **MIM®** technology, combine contemporary design with the comfort for the patient. The anatomical 80 mesh gauge pad fits perfectly on the tooth for an excellent reliability.


Packs of 10

MESIAL FUNNELED ENTRY WAY

FACILITATES THE INSERTION OF THE ARCHWIRE

OCCLUSO-GINGIVAL PARALLEL PLANES

PROVIDE A SECURE GRIP WITH TWEEZERS

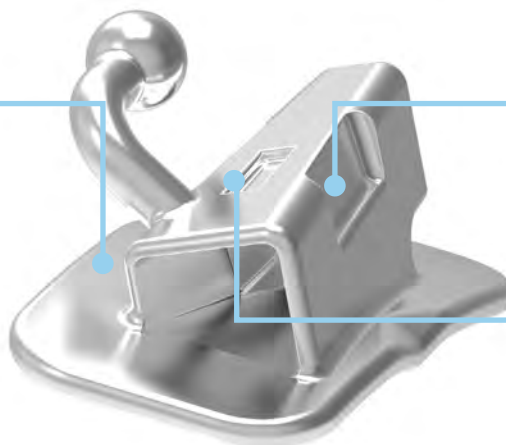

VESTIBULAR INDENT

FACILITATES PROPER POSITION ON THE TOOTH BONDING PHASE

TORQUE-IN-BASE

MINIMIZES OCCLUSAL CONTACTS







IDEAL LENGTH

MAXIMIZES ROTATION CONTROL



ANATOMICAL MESH PAD



WITH INTERCUSPID INDENT, 80 MESH AND FDI IDENTIFICATION

CODE				PRICE			
				torque	rotat.	slot	direct bonding tubes
F8311-16	7 6 6 7	  		-14°	+8°	.022"	L
F8311-26				-14°	+8°	.022"	R
F8311-46	7 6 6 7	  		-20°	+8°	.022"	R
F8311-36				-20°	+8°	.022"	L

**INTRO
PROMO:**
\$24.00

MEMORIA® PREFORMED ARCHWIRES

In a superelastic nickel-titanium base alloy, these archwires provide gentle and continuous forces, reduce treatment time and improve patient comfort. The highly polished surface eliminates any bracket slot friction. Their non-deformability and superelasticity, which are expressed in the continuous force delivery over a wide deflection range, have been tested on thousands of treatments. Available in small, medium and large shape, both upper and lower, they represent the ideal answer to achieve your treatment goals.

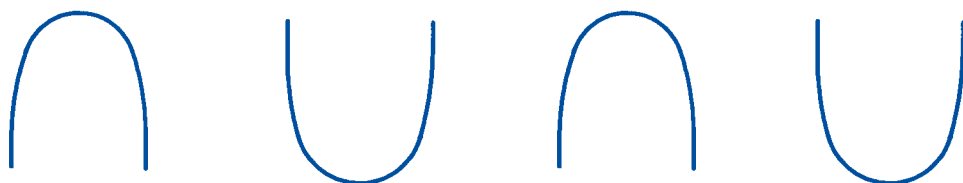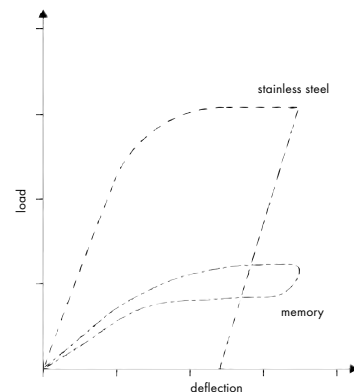


CODE	medium	pkg.	PRICE	CODE	medium	pkg.	PRICE	CODE	medium	pkg.	PRICE	CODE	medium	pkg.	PRICE	Ø
upper				lower				upper				lower				inch
C5910-12	10		\$34.10	C5950-12	10		\$34.10	C5910-012	25		\$80.90	C5950-012	25		\$80.90	.012
C5910-14	10		\$34.10	C5950-14	10		\$34.10	C5910-014	25		\$80.90	C5950-014	25		\$80.90	.014
C5910-16	10		\$34.10	C5950-16	10		\$34.10	C5910-016	25		\$80.90	C5950-016	25		\$80.90	.016
C5910-18	10		\$34.10	C5950-18	10		\$34.10	C5910-018	25		\$80.90	C5950-018	25		\$80.90	.018
C5910-20	10		\$34.10	C5950-20	10		\$34.10	C5910-020	25		\$80.90	C5950-020	25		\$80.90	.020
																inch
C5911-16	10		\$39.40	C5951-16	10		\$39.40	C5911-016	25		\$89.30	C5951-016	25		\$89.30	.016x.016
C5912-16	10		\$39.40	C5952-16	10		\$39.40	C5912-016	25		\$89.30	C5952-016	25		\$89.30	.016x.022
C5912-17	10		\$39.40	C5955-17	10		\$39.40	C5912-017	25		\$89.30	C5955-017	25		\$89.30	.017x.025
C5911-18	10		\$39.40	C5951-18	10		\$39.40	C5911-018	25		\$89.30	C5951-018	25		\$89.30	.018x.018
C5915-18	10		\$39.40	C5955-18	10		\$39.40	C5915-018	25		\$89.30	C5955-018	25		\$89.30	.018x.025
C5915-19	10		\$39.40	C5955-19	10		\$39.40	C5915-019	25		\$89.30	C5955-019	25		\$89.30	.019x.025
C5915-21	10		\$39.40	C5955-21	10		\$39.40	C5915-021	25		\$89.30	C5955-021	25		\$89.30	.021x.025

**SPECIAL
BUY 3
GET 1
FREE**

STAINLESS STEEL PREFORMED ARCHWIRES


Extra spring hard wire, round or rectangular, diamond drawn for leveling of surfaces. Available medium shape, both upper and lower, they represent the ideal answer to achieve your treatment goals.

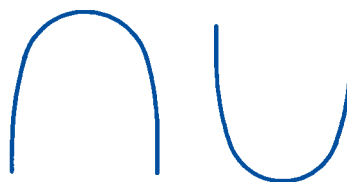

CODE	medium	PRICE	CODE	medium	PRICE	Ø	pkg.
upper			lower			inch	
C3110-14		\$19.40	C3150-14		\$19.40	.014	10
C3110-16		\$19.40	C3150-16		\$19.40	.016	10
C3110-18		\$19.40	C3150-18		\$19.40	.018	10
C3110-20		\$19.40	C3150-20		\$19.40	.020	10
							inch
C3111-16		\$24.90	C3151-16		\$24.90	.016x.016	10
C3112-16		\$24.90	C3152-16		\$24.90	.016x.022	10
C3112-17		\$24.90	C3152-17		\$24.90	.017x.025	10
C3112-18		\$24.90	C3152-18		\$24.90	.018x.025	10
C3112-19		\$24.90	C3152-19		\$24.90	.019x.025	10
C3112-21		\$24.90	C3152-21		\$24.90	.021x.025	10

STRAIGHT ROUND WIRE

An alloy properly designed for forming arches and springs in fixed orthodontic appliances. It is diamond drawn to obtain resistant surfaces. Available in three tempers: hard, spring hard and extra spring hard to be adapted to any therapeutic application.


CODE Hard		pkg.	PRICE
1400/1600 N/mm ²	inch		
C8080-36	.036	20 wires	\$38.90
C8080-40	.040	20 wires	\$38.90
C8080-45	.045	20 wires	\$38.90
C8080-51	.051	20 wires	\$38.90

SPECIAL:
BUY 3 GET
1 FREE!
MIX & MATCH
BY GROUP



STRAIGHT RECTANGULAR WIRE



An alloy properly designed for orthodontic applications. Diamond drawn to obtain perfectly square angles and levelled surfaces. Available in two tempers: spring hard and extra spring hard, it is tested and quality controlled in the laboratory.

CODE spring hard		pkg.	PRICE
1800/2000 N/mm ²	inch		
C8022-16	.016x.022	10 wires	\$44.40
C8025-17	.017x.025	10 wires	\$44.40
C8025-18	.018x.025	10 wires	\$44.40
C8025-19	.019x.025	10 wires	\$44.40


SPECIAL:
BUY 3 GET
1 FREE!
MIX & MATCH
BY GROUP



STRAIGHT TWIST WIRE

CODE			pkg.	PRICE
		inch		
C8580-15		.015	10 wires	\$38.90
C8580-17		.0175	10 wires	\$38.90
C8580-19		.0195	10 wires	\$38.90
C8580-21		.0215	10 wires	\$38.90

SPECIAL:
BUY 3 GET
1 FREE!
MIX & MATCH
BY GROUP

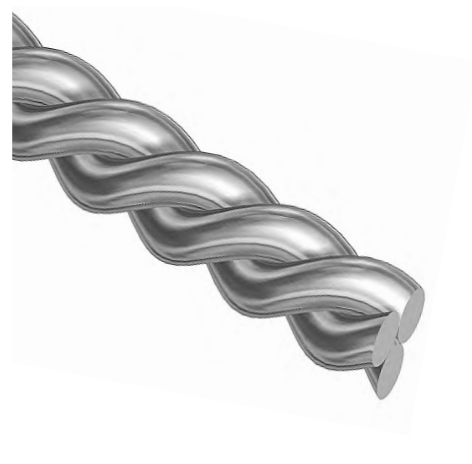

"AUSTRALIAN" WIRE SPOOLS

It is diamond drawn for perfect levelling of surfaces.
The temper is laboratory tested and quality controlled.
Spools, 8 m

CODE extra spring hard	Ø	pkg.	PRICE
2200/2400 N/mm ²	inch		
C0420-12	.012	8 m	\$37.00
C0420-14	.014	8 m	\$37.00
C0420-16	.016	8 m	\$37.00
C0420-18	.018	8 m	\$37.00
C0420-20	.020	8 m	\$37.00



SPECIAL:
BUY 3 GET
1 FREE!
MIX & MATCH
BY GROUP

ELASTOMER PROTECTIVE TUBE

High flexible tube that can easily slip over the arch wire to protect buccal mucosa from damage or irritation. When cut to proper length, it prevents over-closure of spaces.

CODE	Ø	inner	pkg.	PRICE
transparent	mm	inch		
K6940-05	0,5	.018	3 m	\$13.70
K6940-07	0,7	.027	3 m	\$13.70



PALATAL BARS

Palatal bars allow controllable movements of rotation, expansion and contraction of upper first molars. Made of heat-treated wire, 0,9 mm in diameter, to increase their tensile strength, they are easy to remove from the tubes for activation.

ROTATION



EXPANSION



SPECIAL:
BUY 2
5 PACKS
GET ONE
FREE

CODE	mm	pkg.	PRICE
C0570-34	34	5	\$21.50
C0570-36	36	5	\$21.50
C0570-38	38	5	\$21.50
C0570-40	40	5	\$21.50
C0570-42	42	5	\$21.50
C0570-44	44	5	\$21.50
C0570-46	46	5	\$21.50
C0570-48	48	5	\$21.50
C0570-50	50	5	\$21.50
C0570-52	52	5	\$21.50



PEDODONTIC BANDS UNIVERSAL

DESIGNED TO MEET THERAPEUTIC NEEDS
OF PEDIATRIC ORTHODONTICS BY ALLOWING EARLY TREATMENT
OF PATIENTS WITH MIXED OR DECIDUOUS TEETH



SHAPE AND HEIGHT

DESIGNED FOR THE ANATOMY
OF DECIDUOUS MOLARS

EASY FIT

TO THE PYRAMID-SHAPE
OF PRIMARY MOLARS



LASER-MARKING OF SIZES

IT ALLOWS THE STERILIZATION AT
HIGH TEMPERATURES
AND AN OPTIMAL VIEW OF
SIZING OVER TIME

USEFUL IN THE CONSTRUCTION

OF PALATAL EXPANDERS
AND SPACE MAINTAINERS



FREE PREWELDING SERVICE

CUSTOM PREWELDING SERVICE
TO TUBE, BRACKET OR ACCESSORY
WITH ITEM CODE G, CONTACT
LEONEAMERICA FOR ASSISTANCE

PEDODONTIC BANDS UNIVERSAL DECIDUOUS MOLARS UPPER AND LOWER

Available in 10 sizes, either for upper and lower, with LASER etched identification number. The medium temper of biomedical steel used to manufacture the bands allows the adaptation to the pyramid shape of the primary molars by minimizing the problems of detachment. Necessary for the construction of palatal expanders and space maintainers.

**Always specify the size
when ordering e.g.: E6100-05**



CODE		pkg.	PRICE
E6100-00	Upper - PU sizes 1-10	5	\$22.40
E6500-00	Lower - PL sizes 1-10	5	\$22.40

E6165-96 ASSORTMENT PEDODONTIC BANDS UNIVERSAL - KIT

The kit is composed of 5 universal pedodontic bands each upper and lower size, total 100 pcs. The tray is not autoclavable.

PRICE \$175.20



CALIBRA® ANATOMICAL BANDS

MADE OF HARD HEAT-TREATED STAINLESS STEEL

GRADUATED THICKNESS
BY VARIATION IN TOOTH SIZE

**INNER SURFACE
ROUGH FINISH**

FOR MAXIMUM STABILITY
AND INCREASED RETENTION



ANATOMICALLY CORRECT
WIDE RANGE OF SIZES

LASER-MARKING OF SIZES
ALLOWS THE STERILIZATION AT
HIGH TEMPERATURES AND A
DURABLE VIEW OF SIZING



FREE PREWELDING SERVICE

CUSTOM PREWELDING SERVICE TO ANY TUBE, BRACKET OR ACCESSORY WITH ITEM CODE G. READY DELIVERY OF STANDARD PREWELDED BANDS FOR THE MOST COMMONLY USED TECHNIQUES. CONTACT LEONEAMERICA FOR ASSISTANCE.



WIDE CHOICE
OF ASSORTMENTS
IN AUTOCLAVABLE TRAYS

SEE THE CONVERSION CHARTS AT PAGES 65 AND 66 OF **CLINICAL PRODUCTS** SECTION

UPPER AND LOWER BICUSPIDS

Available in 32 sizes either for maxillary or mandibular arch, with LASER etched identification number. Their anatomical shape with contoured edges is designed for use either right or left. Thickness of 0,15mm and optional seating lug. The hard temper and a rough inner surface finish considerably improve cement retention.

Always specify the size when ordering e.g.: E5101-11





**SPECIAL:
BUY 2
10 PKGS
GET 1
FREE**

**SPECIAL:
BUY 2
KITS
GET 1
FREE**

CODE		pkg.	PRICE
E5101-00	Upper with lug - U	10	\$21.50
E5501-00	Lower with lug - L	10	\$21.50
E5100-00	Upper without lug - U	10	\$21.50
E5500-00	Lower without lug - L	10	\$21.50
E5151-00	Upper with lug - U	5	\$14.20
E5551-00	Lower with lug - L	5	\$14.20
E5105-00	Upper without lug - U	5	\$14.20
E5505-00	Lower without lug - L	5	\$14.20

REFERENCE THE CHARTS BELOW
FOR AVERAGE SIZING NEEDS

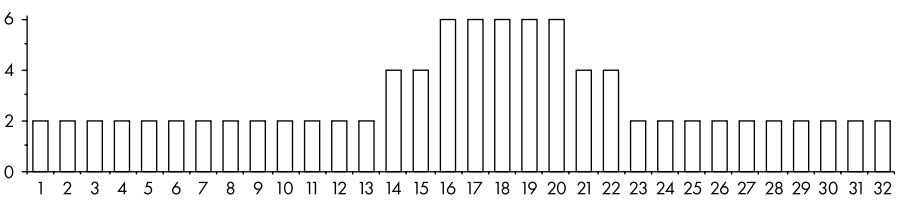ASSORTMENTS	ASSORTED KIT OF 92 BANDS ON AVERAGE USE	CODE	PRICE
quantity each size			
		E5101-92	5 4 4 5 U with lug \$198.70
		E5100-92	5 4 4 5 U without lug \$198.70

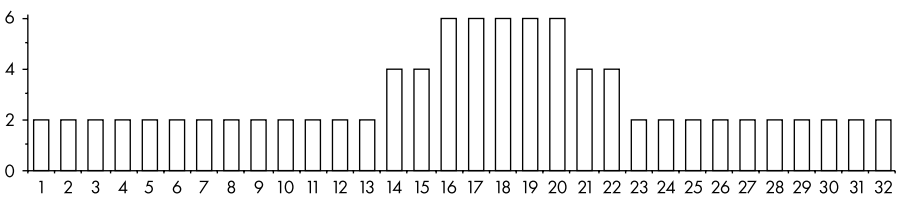ASSORTMENTS	ASSORTED KIT OF 92 BANDS ON AVERAGE USE	CODE	PRICE
quantity each size			
		E5501-92	5 4 4 5 L with lug \$198.70
		E5500-92	5 4 4 5 L without lug \$198.70

UPPER AND LOWER 1ST MOLAR BANDS

Available in 36 right sizes and 36 left sizes either for maxillary or mandibular arch, with LASER etched identification number. The hard temper, the anatomical design, the diversified thickness and the rough inner surface finish make the band selection easy and improve cement retention.

**Always specify the size
when ordering e.g.: E1200-06**



CODE		pkg.	PRICE
E1200-00	Upper right - UR	10	\$21.50
E1300-00	Upper left - UL	10	\$21.50
E1600-00	Lower right - LR	10	\$21.50
E1700-00	Lower left - LL	10	\$21.50
E1205-00	Upper right - UR	5	\$14.20
E1305-00	Upper left - UL	5	\$14.20
E1605-00	Lower right - LR	5	\$14.20
E1705-00	Lower left - LL	5	\$14.20

REFERENCE THE CHARTS BELOW
FOR AVERAGE SIZING NEEDS


K6378-00 RADIO OPAQUE SEPARATORS IN MODULES


Their circular section allows an easy application with special pliers P1141-00 assuring a gradual and comfortable separation of the teeth before the application of the bands. Radio opaque for a secure X-ray visualization.

Pack of 100 modules, 10 rings each

PRICE \$80.00

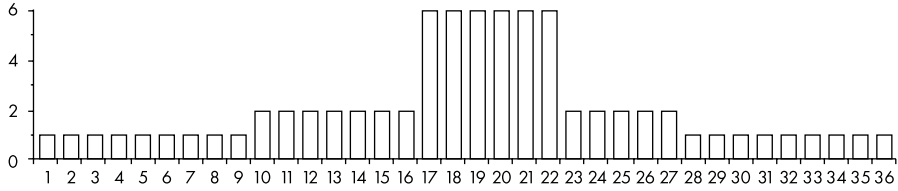

ASSORTMENTS	ASSORTED KIT OF 78 BANDS ON AVERAGE USE	CODE	PRICE
quantity each size			
		E1200-92	<u>6</u> UR R \$168.50
		E1300-92	<u>6</u> UL L \$168.50
		E1600-92	<u>6</u> LR R \$168.50
		E1700-92	<u>6</u> LL L \$168.50

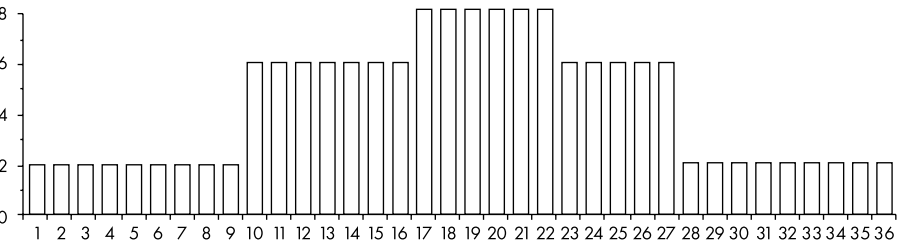ASSORTMENTS	ASSORTED KIT OF 156 BANDS ON AVERAGE USE	CODE	PRICE
quantity each size			
		E1200-94	<u>6</u> UR R \$336.80
		E1300-94	<u>6</u> UL L \$336.80
		E1600-94	<u>6</u> LR R \$336.80
		E1700-94	<u>6</u> LL L \$336.80

PREWELDED BANDS

REFERENCE THE CHARTS ON THE PREVIOUS PAGE FOR AVERAGE SIZING NEEDS



EDGEWISE TECHNIQUE - UPPER

Tube .022"x.028".

Round tube .051".

(G8025-02 - G8025-03)

CODE		pkg.	PRICE
E1420-00	Right UR	10	\$56.20
E1430-00	Left UL	10	\$56.20
E1425-00	Right UR	5	\$31.00
E1435-00	Left UL	5	\$31.00



EDGEWISE TECHNIQUE - LOWER

Tube .022"x.028".

Round tube .051".

(G8025-03 - G8025-02)

CODE		pkg.	PRICE
E1460-00	Right LR	10	\$56.20
E1470-00	Left LL	10	\$56.20
E1465-00	Right LR	5	\$31.00
E1475-00	Left LL	5	\$31.00

ANDREWS ROTH BOSTON ALEXANDER TECHNIQUE - UPPER

Tube .022"x.028".

Torque -10°, rotation +8°.

Round tube .051".

(G8325-02 - G8325-03)



CODE		pkg.	PRICE
E1520-00	Right UR	10	\$56.20
E1530-00	Left UL	10	\$56.20
E1525-00	Right UR	5	\$31.00
E1535-00	Left UL	5	\$31.00

ANDREWS ROTH BOSTON ALEXANDER TECHNIQUE - LOWER

Tube .022"x.028".

Torque -25°, rotation +8°.

Round tube .051".

(G8325-06 - G8325-07)






CODE		pkg.	PRICE
E1560-00	Right LR	10	\$56.20
E1570-00	Left LL	10	\$56.20
E1565-00	Right LR	5	\$31.00
E1575-00	Left LL	5	\$31.00




CONVERTIBLE MIM® BUCCAL TUBES

The laser welded convertible cap prevents from any detachment and it is easy removable with the special Leone instrument P1090-00 (Page 55). Every tube is available either in the DB version, or prewelded to CALIBRA®.

Packs of 10. Packs of 100 on request.

CODE		lo		PRICE				
			torque	rect. tube rot.	diam. tube rond	rect. tube slot	weldable tubes	
 auxiliary tube .018"x.025"	G8184-22		0°	0°	.045"occl.	.018"	R	\$35.70
	G8184-23	6/6	0°	0°	.045"occl.	.018"	L	\$35.70
	G8124-22	6/6	0°	0°	.045"occl.	.022"	R	\$35.70
	G8124-23		0°	0°	.045"occl.	.022"	L	\$35.70
	G8185-12		0°	0°	.051"occl.	.018"	R	\$35.70
	G8185-13	6/6	0°	0°	.051"occl.	.018"	L	\$35.70
	G8125-12	6/6	0°	0°	.051"occl.	.022"	R	\$35.70
	G8125-13		0°	0°	.051"occl.	.022"	L	\$35.70
	G8484-22		-10°	+8°	.045"occl.	.018"	R	\$35.70
	G8484-23	6/6	-10°	+8°	.045"occl.	.018"	L	\$35.70
	G8424-22		-10°	+8°	.045"occl.	.022"	R	\$35.70
	G8424-23		-10°	+8°	.045"occl.	.022"	L	\$35.70
	G8485-12		-10°	+8°	.051"occl.	.018"	R	\$35.70
	G8485-13	6/6	-10°	+8°	.051"occl.	.018"	L	\$35.70
	G8425-12		-10°	+8°	.051"occl.	.022"	R	\$35.70
	G8425-13		-10°	+8°	.051"occl.	.022"	L	\$35.70
 auxiliary tube .018"x.025"	G8170-16		0°	0°		.018"	R	\$35.70
	G8170-17	6/6	0°	0°		.018"	L	\$35.70
	G8110-16	6/6	0°	0°		.022"	R	\$35.70
	G8110-17		0°	0°		.022"	L	\$35.70
	G8470-02		-14°	+8°		.018"	R	\$35.70
	G8470-03	6/6	-14°	+8°		.018"	L	\$35.70
	G8410-02		-14°	+8°		.022"	R	\$35.70
	G8410-03		-14°	+8°		.022"	L	\$35.70
	G8470-06		-25°	+8°		.018"	R	\$35.70
	G8470-07	6/6	-25°	+8°		.018"	L	\$35.70
	G8410-06		-25°	+8°		.022"	R	\$35.70
	G8410-07		-25°	+8°		.022"	L	\$35.70
 auxiliary tube .018"x.025"	G8184-12		0°	0°	.045"occl.	.018"	R	\$35.70
	G8184-13	6/6	0°	0°	.045"occl.	.018"	L	\$35.70
	G8124-12	6/6	0°	0°	.045"occl.	.022"	R	\$35.70
	G8124-13		0°	0°	.045"occl.	.022"	L	\$35.70
	G8185-22		0°	0°	.051"occl.	.018"	R	\$35.70
	G8185-23	6/6	0°	0°	.051"occl.	.018"	L	\$35.70
	G8125-22	6/6	0°	0°	.051"occl.	.022"	R	\$35.70
	G8125-23		0°	0°	.051"occl.	.022"	L	\$35.70
	G8484-12		-10°	+8°	.045"occl.	.018"	R	\$35.70
	G8484-13	6/6	-10°	+8°	.045"occl.	.018"	L	\$35.70
	G8424-12		-10°	+8°	.045"occl.	.022"	R	\$35.70
	G8424-13		-10°	+8°	.045"occl.	.022"	L	\$35.70
	G8485-22		-10°	+8°	.051"occl.	.018"	R	\$35.70
	G8485-23	6/6	-10°	+8°	.051"occl.	.018"	L	\$35.70
	G8425-22		-10°	+8°	.051"occl.	.022"	R	\$35.70
	G8425-23		-10°	+8°	.051"occl.	.022"	L	\$35.70
	G8424-32	6/6	-14°	+10°	.045"occl.	.022"	R	\$35.70
	G8424-33		-14°	+10°	.045"occl.	.022"	L	\$35.70

CONVERTIBLE MIM® BUCCAL TUBES

CODE							PRICE	
			torque	rect. tube rot.	diam. tube rond	rect. tube slot	weldable tubes	
	G8184-02		0°	0°	.045"ging.	.018"	R	\$35.70
	G8184-03	$\frac{6}{6}$	0°	0°	.045"ging.	.018"	L	\$35.70
	G8124-02		0°	0°	.045"ging.	.022"	R	\$35.70
	G8124-03		0°	0°	.045"ging.	.022"	L	\$35.70
	G8185-02		0°	0°	.051"ging.	.018"	R	\$35.70
	G8185-03	$\frac{6}{6}$	0°	0°	.051"ging.	.018"	L	\$35.70
	G8125-02		0°	0°	.051"ging.	.022"	R	\$35.70
	G8125-03		0°	0°	.051"ging.	.022"	L	\$35.70
	G8484-02		-10°	+8°	.045"ging.	.018"	R	\$35.70
	G8484-03	$\frac{6}{6}$	-10°	+8°	.045"ging..	.018"	L	\$35.70
	G8424-02		-10°	+8°	.045"ging.	.022"	R	\$35.70
	G8424-03		-10°	+8°	.045"ging.	.022"	L	\$35.70
	G8485-02		-10°	+8°	.051"ging.	.018"	R	\$35.70
	G8485-03	$\frac{6}{6}$	-10°	+8°	.051"ging.	.018"	L	\$35.70
	G8425-02		-10°	+8°	.051"ging.	.022"	R	\$35.70
	G8425-03		-10°	+8°	.051"ging.	.022"	L	\$35.70
	G8424-16		-20°	+8°	.045"ging.	.022"	R	\$35.70
	G8424-17		-20°	+8°	.045"ging.	.022"	L	\$35.70
	G8484-06	$\frac{6}{6}$	-25°	+8°	.045"ging.	.018"	R	\$35.70
	G8484-07		-25°	+8°	.045"ging.	.018"	L	\$35.70
	G8424-06		-25°	+8°	.045"ging.	.022"	R	\$35.70
	G8424-07		-25°	+8°	.045"ging.	.022"	L	\$35.70
	G8485-06		-25°	+8°	.051"ging.	.018"	R	\$35.70
	G8485-07	$\frac{6}{6}$	-25°	+8°	.051"ging.	.018"	L	\$35.70
	G8425-06		-25°	+8°	.051"ging.	.022"	R	\$35.70
	G8425-07		-25°	+8°	.051"ging.	.022"	L	\$35.70
	G8181-02		0°	0°		.018"	R	\$35.70
	G8181-03	$\frac{6}{6}$	0°	0°		.018"	L	\$35.70
	G8121-02		0°	0°		.022"	R	\$35.70
	G8121-03		0°	0°		.022"	L	\$35.70
	G8481-02		-10°	+8°		.018"	R	\$35.70
	G8481-03	$\frac{6}{6}$	-10°	+8°		.018"	L	\$35.70
	G8421-02		-10°	+8°		.022"	R	\$35.70
	G8421-03		-10°	+8°		.022"	L	\$35.70
	G8421-16		-20°	+8°		.022"	R	\$35.70
	G8421-17	$\frac{6}{6}$	-20°	+8°		.022"	L	\$35.70
	G8481-06		-25°	+8°		.018"	R	\$35.70
	G8481-07	$\frac{6}{6}$	-25°	+8°		.018"	L	\$35.70
	G8421-06		-25°	+8°		.022"	R	\$35.70
	G8421-07		-25°	+8°		.022"	L	\$35.70



MIM® LINGUAL TUBE WITH BALL HOOK

- G2846-02** Right
G2846-03 Left
G2847-02 Right, rotation -8°
G2847-03 Left, rotation -8°



The curved flanges adapt to contours of both first and second molars. For use with double wires: .036".

Packs of 10 **PRICE \$31.20**

SPECIAL:
 BUY 2 PACKS &
 PAY \$24.00
 EACH MIX &
 MATCH



LINGUAL BUTTONS

G2860-00
 Flat for centering on molars



G2865-00
 Curved for mesial placement on molars



G2870-00
 Deep curved for bicuspid

Packs of 10 **PRICE \$26.30**



F2810-00
 Direct bond eyelet

Packs of 10 **PRICE \$26.00**

SPECIAL:
 BUY 6 PACKS,
 PAY \$19.00/PK
 MIX & MATCH



LINGUAL CLEATS LOW PROFILE

- G2824-00** For molars
G2822-00 For cuspids and bicuspid



Packs of 100 **PRICE \$138.30**

SPECIAL:
 BUY 2 AT
\$94.00 EACH



CRIMPABLE ARCH STOP

It can be inserted without removing the archwire. Once placed into the mouth, crimp it with a heavy cutter or universal pliers. Suggested for use with coil springs too.

- G9225-18** arches .018"x.025"
G9228-22 arches .022"x.028"

Packs of 10 **PRICE \$21.80**



MIM® SURGICAL HOOKS

For intermaxillary post-surgical fixation and tractions. Round edges. Easy to place with plier P1911-00 (Page 55).

CODE	CODE	←→	pkg.	PRICE
for arch .018"	for arch .022"	mm		
G9620-18	G9620-22	2,0	10	\$34.10
G9650-18	G9650-22	2,0 open	10	\$34.10

LIGATURE WIRE

Made of stainless steel.



CODE	Ø		pkg.	PRICE
	mm	inch		
K4050-23	0,23	.009	spools of 20 g approx.	\$16.80
K4050-25	0,25	.010		\$16.80
K4050-28	0,28	.011		\$16.80
K4050-30	0,30	.012		\$16.80
K4051-23	0,23	.009	reels of 500 g approx.	\$75.10
K4051-25	0,25	.010		\$75.10
K4051-28	0,28	.011		\$75.10
K4051-30	0,30	.012		\$75.10
K4051-36	0,36	.014		\$75.10

PREFORMED WIRE KOBAYASHI

Once ligature is ready, a traction hook may be formed with the looped soldered end, suitable for any kind of bracket.

CODE	Ø		pkg.	PRICE
	mm	inch		
K4610-30	0,30	.012	100	\$23.90
K4610-36	0,36	.014	100	\$23.90



PREFORMED LIGATURE WIRE

CODE	Ø		pkg.	PRICE
	mm	inch		
K4210-23	0,23	.009	1000	\$38.90
K4210-25	0,25	.010	1000	\$38.90
K4210-28	0,28	.011	1000	\$38.90
K4210-30	0,30	.012	1000	\$38.90



ELASTIC WIRE FOR LIGATURES AND TRACTION

Made of highly elastic, atoxic polyurethane. The rough surface enhances the tightness of ligatures and traction loops.

CODE	Ø	pkg.	PRICE
transparent	mm		
K6125-06	0,6	8 m	\$24.90
K6125-08	0,8	8 m	\$24.90



ELASTIC CHAINS FOR LIGATURES

Manufactured out from high elastic material with a precision injection moulding process. The perfect ring section facilitates the ligature of the brackets. The use of the mini modules allows for a better control of infections.







CODE	CODE	Ø	pkg.	PRICE
grey	transparent	inner		
K6460-13	K6420-13	long	1,30 mm	spools
K6465-13	K6425-13	short	1,30 mm	of 5 m
K6466-13	K6426-13	extra short	1,30 mm	(15 feet)
				\$29.90
				\$29.90
				\$29.90

INTRAORAL ELASTICS

Cut from pure surgical latex tubing. All packages indicate the diameter size and the elastic pull which comes into force when the elastic is stretched three times its diameter. On each bag a colored image is printed, for an immediate identification from young patients.

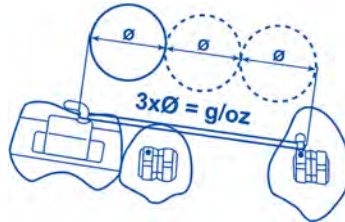Packs of 10 bags, 100 elastics approx. each bag



CODE	CODE	CODE	CODE	CODE	Ø inner		PRICE
brown 1/8" - 3,2 mm	green 3/16" - 4,8 mm	red 1/4" - 6,4 mm	blue 5/16" - 7,9 mm	orange 3/8" - 9,5 mm			
K0901-02	K0902-02	K0903-02	K0904-02	K0905-02	light 2,5 oz - 70g		\$19.20
K0906-04	K0907-04	K0908-04	K0909-04	K0910-04	medium 4,5 oz - 126 g		\$19.20
K0915-06	K0911-06	K0912-06	K0913-06	K0914-06	heavy 6,5 oz - 182 g		\$19.20

ELASTIC LIGATURES

Manufactured out from high elastic material with a precision injection moulding process. The perfect ring section facilitates the ligature of the brackets. The use of the mini modules allows for a better control of infections.



ELASTIC COLORED LIGATURE KIT MINI MODULES

K7351-93 inner diam. 1,0 mm - **small** **PRICE \$87.40**

K7352-93 inner diam. 1,3 mm - **medium** **PRICE \$87.40**

K7353-93 inner diam. 1,6 mm - **large** **PRICE \$87.40**

Packs of 840 rings, 6 mini modules of 10 rings each in 7 colors

K7320-93 ELASTIC LIGATURE KIT SILVER AND TRANSPARENT MINI MODULES


Pack of 840 rings:
84 assorted mini modules of 10 rings each:
24 small, 36 medium, 24 large

PRICE \$87.40



ELASTIC LIGATURE REFILLS

CODE	CODE	CODE	colour	pcs.	PRICE
small inner. diam. 1,0 mm	medium inner. diam. 1,3 mm	large inner. diam. 1,6 mm			
K7351-10A	K7352-13A	K7353-16A	blue	200	\$23.60
K7351-10S	K7352-13S	K7353-16S	pink	200	\$23.60
K7351-10V	K7352-13V	K7353-16V	green	200	\$23.60
K7351-10G	K7352-13G	K7353-16G	yellow	200	\$23.60
K7351-10R	K7352-13R	K7353-16R	red	200	\$23.60
K7351-10E	K7352-13E	K7353-16E	orange	200	\$23.60
K7351-10L	K7352-13L	K7353-16L	purple	200	\$23.60
K7321-10	K7322-13	K7323-16	transpa- rent	200	\$23.60
K7361-10	K7362-13	K7363-16	silver	200	\$23.60

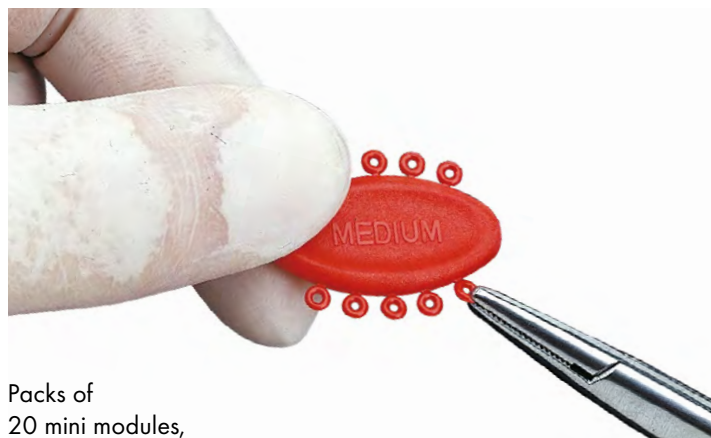

Packs of
20 mini modules,
10 rings each

BONDED RETAINERS

Tooth shaped pads allow the maximum retention strength. Small thickness for a perfect adjustment to the lingual surface, with ease of hygiene and lower patient discomfort.

UPPER RETAINER CENTRAL TO CENTRAL

F3812-93

UPPER RETAINER CENTRAL TO CENTRAL KIT



Assorted kit of 20

PRICE \$121.80

CODE	pkg.		PRICE
	mm		
F3812-04	4	10	\$60.90
F3812-06	6	10	\$60.90
F3812-09	9	10	\$60.90



UPPER RETAINER LATERAL TO LATERAL

F3822-93

UPPER RETAINER LATERAL TO LATERAL KIT

Assorted kit of 20

PRICE \$121.80



CODE	pkg.		PRICE
	mm		
F3822-14	14	10	\$60.90
F3822-16	16	10	\$60.90
F3822-18	18	10	\$60.90
F3822-20	20	10	\$60.90



LOWER RETAINER CUSPID TO CUSPID

F3836-93

LOWER RETAINER CUSPID TO CUSPID KIT

Assorted kit of 10

PRICE \$72.20



CODE	pkg.		PRICE
	mm		
F3836-29	29	2	\$14.40
F3836-31	31	2	\$14.40
F3836-33	33	2	\$14.40
F3836-35	35	2	\$14.40
F3836-37	37	2	\$14.40




SELF-DRILLING ORTHODONTIC MINI IMPLANTS

LeoneAmerica is proud to announce that it is the sole distributor of PSM medical solutions products in United States. In collaboration with leading specialists, PSM has developed the unique Benefit TAD System for temporary fixation in the palate that ensures a high reliable palatal anchorage with Leone RPE4Tads AND Leone CAD/CAM expanders hybrid or completely bone borne. PSM offers high precision titanium medical grade implant screw and a large variability of prefabricated connecting elements and instruments. Full catalog available on request.

**NEW
LINE**

psm
MEDICAL SOLUTIONS

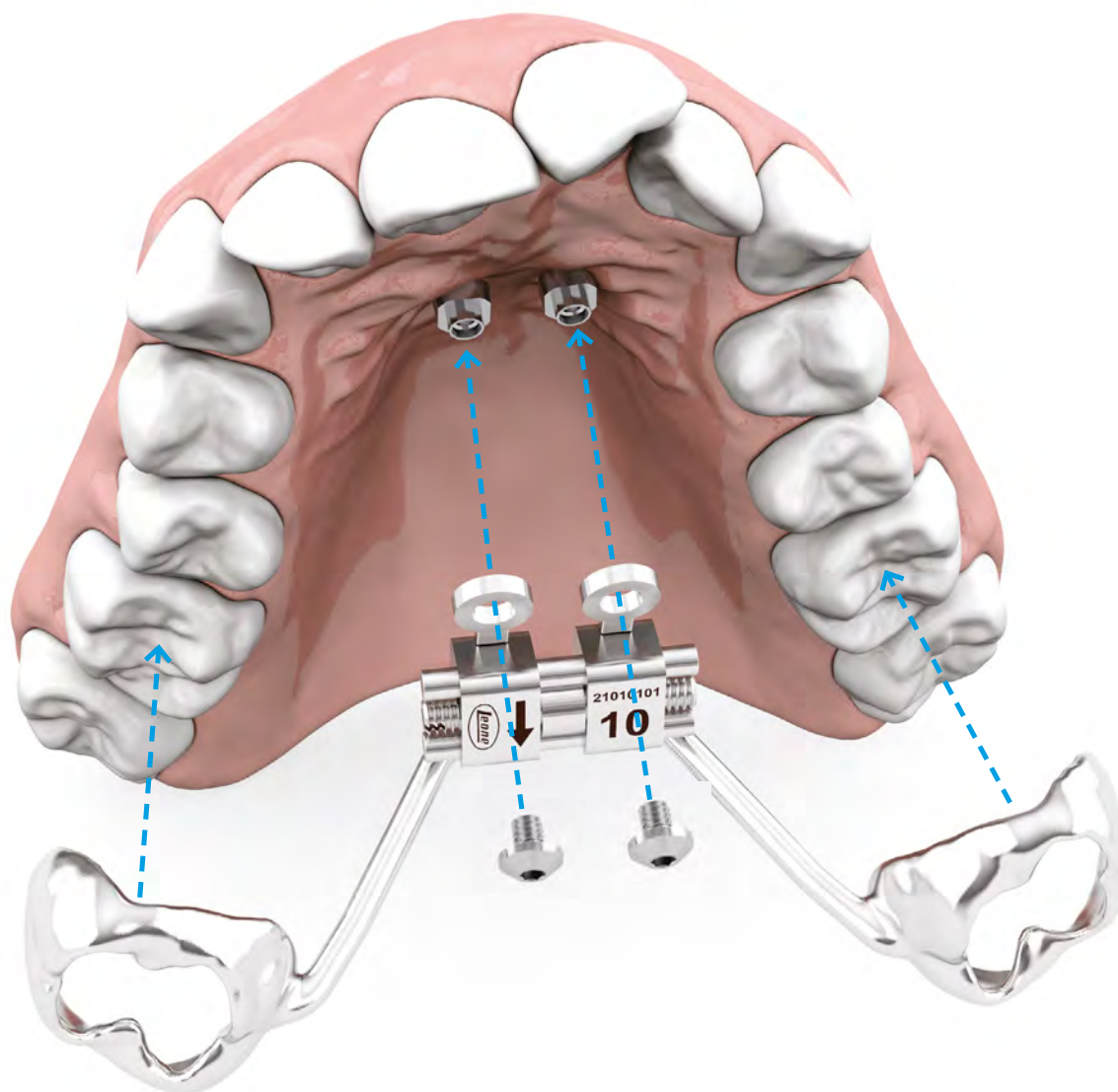

THE TRADEMARKS SHOWN FROM
PAGE 38 TO PAGE 43 ARE NOT
OF LEONE S.P.A. NOR OF LEONE
AMERICA DENTAL PRODUCTS INC.



Call for a PSM pdf
copy of the catalog
or click here

Double hexagon to prevent rotation of abutments.

Internal thread for BENEplate fixation screw or the abutment screw.

BENEFIT

Increased thread depth and denser pitch design ensure maximum anchorage even in cancellous bone.

Continuously sharp thread flanks cut into the bone.

Self-drilling screw tip.

Optimized thread design leads to easy insertion without force (cork screw principal).

Non removable internal abutment screws.




single sterile packaging




Standard
Abutment



Abutment
with loop acc.
to
Dr. S. Renger,
St. Louis,
France

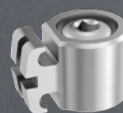

Abutment with
1 or 2 brackets

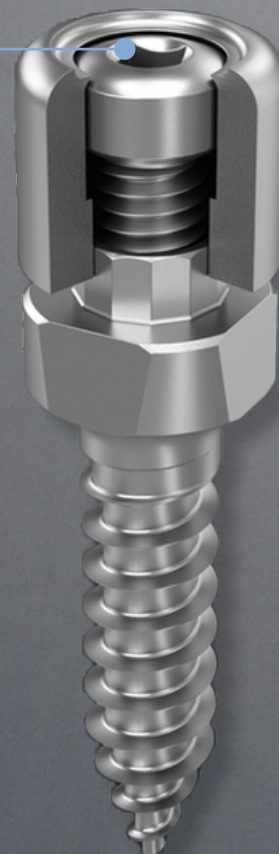

for more detailed informations
download our new BENEFIT catalogue



Skeletal anchorage in Orthodontics

An important requirement for successful orthodontic treatment is stable anchorage. Conventional orthodontic anchorage can be limited by periodontal disease and/or missing teeth. Skeletal anchorage has been proven to be very effective. It can be an essential part of non-compliance therapy and broadens the options available to the orthodontic practitioner.

Mini-Implants are especially advantageous because of the minimally invasive insertion protocol and their relatively low cost.

However, limitations can occur when conventional Mini-Implant Systems are used for orthodontic anchorage. In collaboration with leading specialists, **psm** has developed the unique **BENEFIT®** TAD System for temporary fixation in the palate. **BENEFIT®** offers high precision titanium implant screws and a large variety of pre-fabricated connecting parts such as **BENEplates**, abutments and palatal expanders.


With the introduction of the **BENEFIT®** System **psm** has considerably enhanced the therapeutic diversity of the treatments available for Mesialization, Distalization, RPE as well as uprighting or anchorage of the dentition.

With **BENEFIT®**, **psm** was the first company to offer a full range of skeletal anchorage appliances for the palate – a trend that has become state-of-the-art in Orthodontics.





Distalization with *BENEplates*

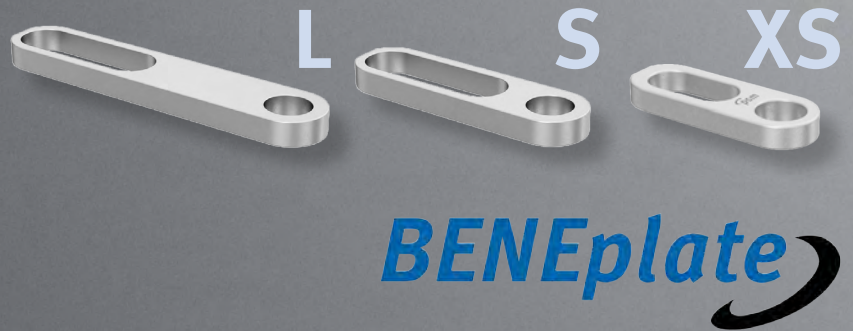


Mesialization with *BENEplates*



Hybrid Hyrax



Starter Kits *BENEFit*®

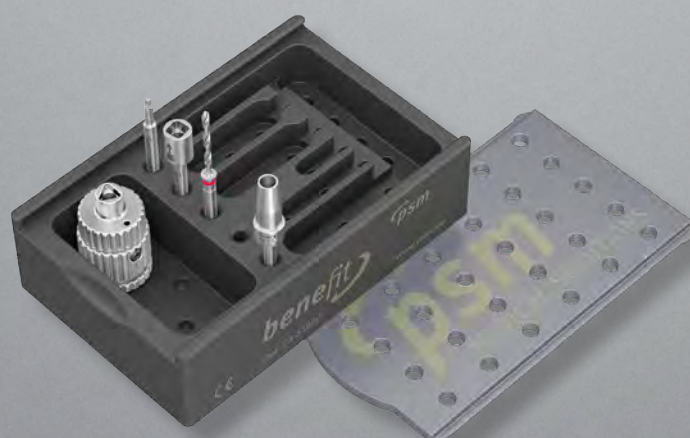

STARTER KIT – M

Examples



Distalization



Mouse Trap



Front Anchorage

For more information check our website www.psm-na.us



Appliance first!

BENEplate DIRECT with prethreaded screw holes guarantee a perfect connection with the **BENEFIT®** DIRECT Screws even in case of a non-axial insertion (up to 15°)

smartlock
Available from Spring 2022

Unthreaded area for atraumatic anchorage within the soft tissue

TX - self holding screw head with increased surface contact

Increased thread depth and denser pitch design ensure maximum anchorage even in cancellous bone

Reduced insertion torque thanks to the specially designed tip and shaft

Precision manufactured thread profile for sharpness and increased self tapping properties

bmX DIRECT




▶ VIDEO



quattro

The **quattro** TAD system offers different screwhead designs which are perfectly adapted to the Orthodontist's appliances used daily. All **quattro** screws feature the identical and well-proven thread design for easy insertion and most secure and satisfactory treatment results.

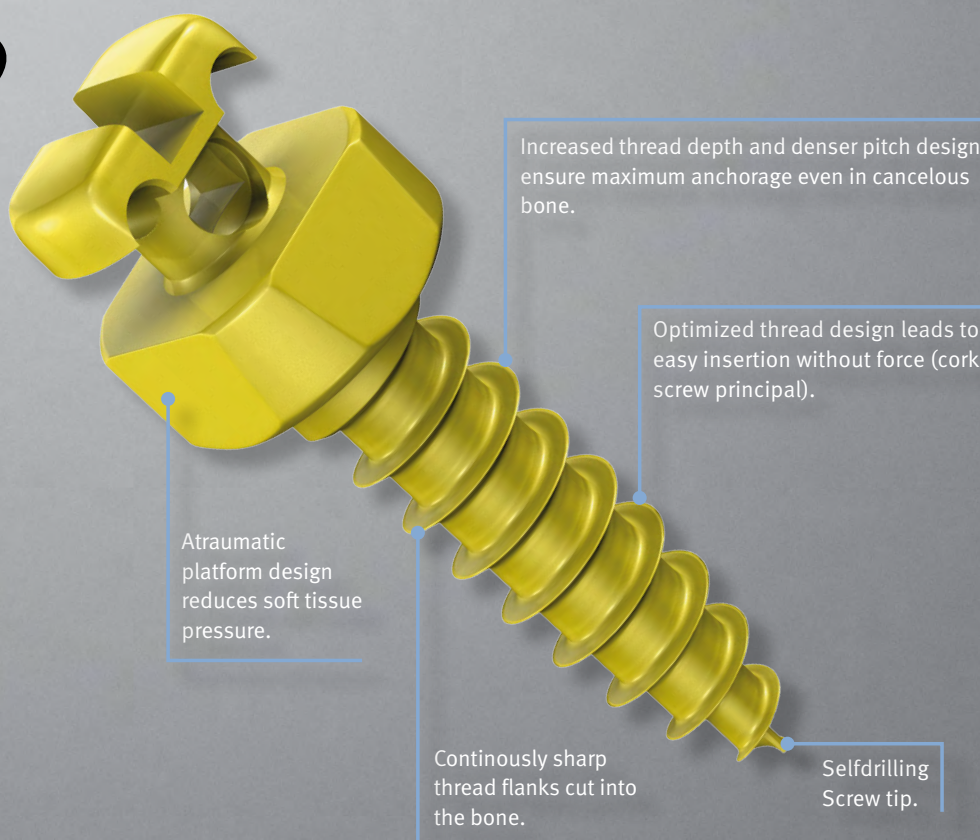

quattro plus



quattro hook



quattro plus V

for more detailed informations
download our quattro catalogue



ORTHODONTIC MINI IMPLANTS

SELF-DRILLING

The special design of the devices makes the insertion possible without the need of any pre drilling to prepare the implant site. The endosseous portion of the self-drilling mini implant features one apical portion (2 mm long) drill shaped in order to punch the bone and one threaded portion (1,75 mm in diameter and 6 mm long) with self-tapping coil for anchorage to the bone.

CLOSED DUAL SLOT





NO OSSEOINTEGRATION

MADE OF SURGICAL
STAINLESS STEEL, EASY
REMOVABLE

low head high head



**NON-STERILE PACKAGE
WITH IMPLANT CARRIER
AND VIAL CAP**
COLOR CODED

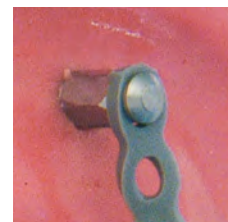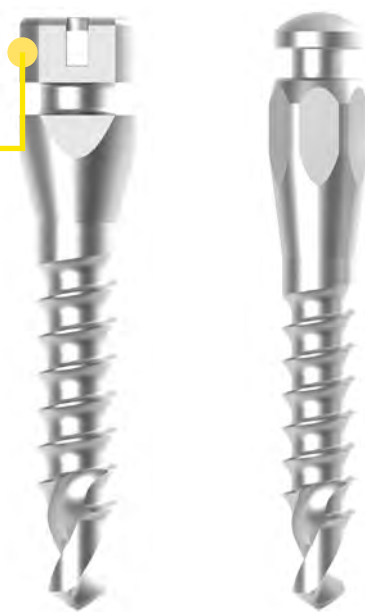

GROOVE
FOR EASY APPLICATION
OF ELASTICS
AND ELASTIC CHAINS

SELF-DRILLING TIP
DRILL SHAPED

INFORMATION FOR DISTRIBUTORS OF ORTHODONTIC MINI IMPLANTS, "INTENDED USE", "RESPONSIBILITY", "SURVEILLANCE"



The 93/42EEC directive on medical devices is the official reference that dictates the regulations for marketing medical devices. The directive provides indications for all the phases of existence for the device (from the project phase through the traceability system, and surveillance), and it identifies all the characters who have to comply with the directive itself, which includes not only the manufactures, but also the distributors, the buyers, and even the users. As for the responsibilities of the single competence, Leone S.p.A. recommends to its direct clients, dental depots and exclusive dealers to follow and maintain the indications, warnings, and information for the univocal identification of the medical devices, as provided by the manufacturer on the labels, during all the marketing phases. With specific regard to Class IIB implantable products, all dental depots and exclusive dealers of Leone S.p.A. are required to keep records of the distribution of medical devices as of traceability available for verification, in case of need to trace back a product or its user in a univocal way.

SELF-DRILLING ORTHODONTIC MINI IMPLANTS WITH CLOSED DUAL SLOT

The head of the mini implant is available in two transmucosal heights and has a **closed dual slot** for the application of wire and/or spring segments. The **closed dual slot** section is .018"x.025" and can accommodate round or rectangular wires. A minimal invasive surgical procedure is needed: in one phase, transmucosal with very good primary stability for immediate application of the orthodontic traction.

Mini implants are supplied non-sterile (for the sterilization process see the instruction for use in the package), please use with the proper adapter for handpiece or with the screwdriver for hand insertion.

Packs of 1

CODE		Ø		PRICE
		mm	mm	
002-1708P03	low head			
		1,75	8	\$80.90
002-1708P04	high head			
		1,75	8	\$80.90

**CLOSED
DUAL
SLOT**



080-1002-01

HANDPIECE ADAPTER FOR MINI IMPLANTS WITH CLOSED DUAL SLOT

For use with the handpiece for the insertion of self-drilling mini implants with closed dual slot also in difficult access areas.

PRICE \$68.30



080-1000-02

SCREWDRIVER FOR MINI IMPLANTS WITH CLOSED DUAL SLOT

PRICE \$169.60

080-1000-23

REFILL TIP



PRICE \$48.30




SELF-DRILLING ORTHODONTIC MINI IMPLANTS


The head of the self-drilling mini implant is available in two transmucosal heights, it is mushroom shaped to facilitate the application of elastic ligatures and coil springs. One stage, transmucosal and minimally invasive surgical technique is needed, with optimal primary stability for the immediate load of the orthodontic appliance. Self-drilling mini implants are supplied not-sterile (for the sterilization process see the instruction for use in the package), to use with the special handpiece adapter. The use of a surgical micromotor is recommended.

Packs of 1

CODE		Ø		PRICE
		mm	mm	
002-1708P01	low head 	1,75	8	\$80.90
002-1708P02	high head 	1,75	8	\$80.90

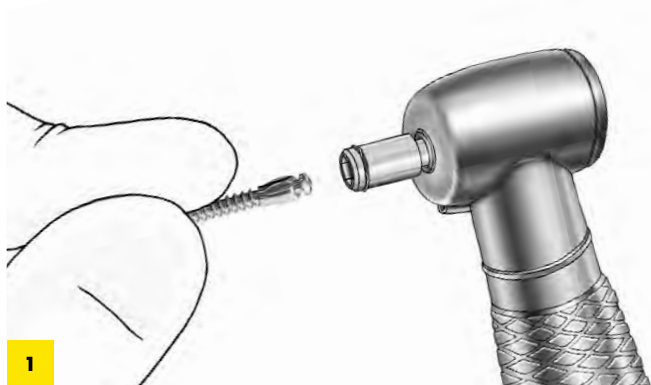

Once sterilized insert the head of the mini implant to the handpiece adapter Cat. 080-1002-00 (for Self-Drilling Orthodontic Mini implant) or Cat. 080-1002-01 (for Self-Drilling mini implants with closed dual slot) connected to a surgical micromotor.

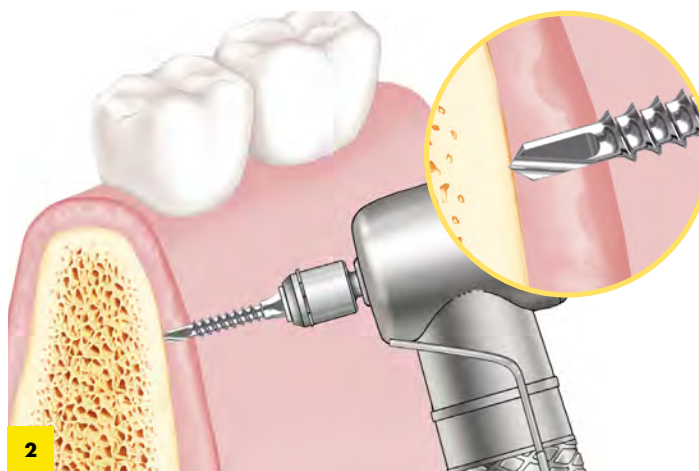

Lean the drilling tip of the mini implant on the attached mucosa in the selected area. Exert a pressure gentle enough to punch the soft tissues till the cortical bone is reached.



Set the maximum speed on micromotor to 300 rpm and torque value to 40 Ncm. Start drilling the mini implant while irrigating abundantly with water. Immediately after having reached the cortical bone and as soon as the presence of a softer bone tissue has been detected, reduce the micromotor maximum speed to 50 rpm. Complete the insertion under these parameters.



080-1002-00 HANDPIECE ADAPTER

For use with the contra-angle for the insertion of mini implants also in difficult access areas.

PRICE \$43.60

MINI IMPLANT EXTRACTION.

AT THE END OF TREATMENT, AND NEVER LATER THAN 6 MONTHS, THE MINI IMPLANTS HAVE TO BE REMOVED. TO REMOVE MINI IMPLANTS, THE SAME INSTRUMENTS ALREADY EMPLOYED FOR THE INSERTION CAN BE USED BY TURNING THEM COUNTERCLOCKWISE.



The pictures and information described in this page have informative purpose only and they are not intended to furnish any therapeutic advice for the dental surgeon, the dentist, the dental technician or the patient. Leone Spa expressly disclaims any liability neither will furnish any warranty regarding the exactness or the conformance of the information furnished in these pages.

FACE MASKS






Used for treatment of maxillary insufficiencies, prognathism of the mandible and cleft palate. The reverse orthopaedic pull between maxilla and mandible facilitates the subsequent orthodontic treatment. Made of stainless steel wire and non toxic acrylic, it can be adapted to perform any individual size adjustment, either vertical or horizontal. All attachments are locked by set screws in metal housings. Face masks with zygomatic holds feature an adjustable strap band.

Packs of: 1 face mask, 1 screw wrench, instructions for use in a sleeve for patient



REPLACEMENT PARTS

M0777-00	Chin pad	Pack of 1	PRICE \$15.00
M0778-00	Forehead pad	Pack of 1	PRICE \$15.00
M0779-00	Zygomatic hold pad	Pack of 2	PRICE \$15.00

	CODE	CODE	CODE	CODE		PRICE
WITH HORIZONTAL ADJUSTMENT	BEIGE	RED	YELLOW	GREEN		
	M0774-00	M0774-00R	M0774-00G	M0774-00V	SMALL	\$68.30
	M0776-00	M0776-00R	M0776-00G	M0776-00V	LARGE	\$68.30
	M0773-00				MINI	\$49.10
UNIVERSAL ADJUSTABLE						
	M0775-00	M0775-00R	M0775-00G	M0775-00V		\$81.90
WITH VERTICAL ADJUSTMENT						
	M0774-01				SMALL	\$68.30
	M0776-01				LARGE	\$68.30
WITH ZYGOMATIC HOLDS						
	M0771-00					\$62.00
DYNAMIC FACE MASK						
VERTICAL ADJUSTMENT	BLUE	RED				
	M0772-01	M0772-01R				\$112.40

ok ORTHODONTICS
FOR KIDS

SPECIAL:
BUY 2
MIX & MATCH
GET 15% OFF

CODE =inner diameter	CODE =inner diameter			PRICE
3/8" - 9,5 mm	1/2" - 12,7 mm			
K0935-08	K0936-08	strong	8 oz - 224 gr	\$4.20



HEADGEAR FOR SAFETY MODULES

High pull. Use with safety release spring modules only. It fits onto safety hooks of spring modules.

Made of hypoallergenic and washable material with comfortable padding. Available in 3 sizes and 2 colors to provide patient's psychological compliance.

Packs of 1

CODE		PRICE
blue		
M0805-00	small	\$21.50
M0805-10	medium	\$21.50
M0805-20	large	\$21.50
pink		
M0805-00S	small	\$21.50
M0805-10S	medium	\$21.50
M0805-20S	large	\$21.50



NECK PAD FOR SAFETY MODULES

Low pull. Same as headgear both for use with safety modules and materials.

Available in 2 colors to provide patient's psychological compliance.

Packs of 10

M0800-00S Pink

PRICE \$62.50

M0800-00 Blue

PRICE \$62.50






EXTRAORAL ELASTICS

Cut from pure surgical latex.

Pack of 1 bag

with 100 elastics approx.



CODE		CODE		CODE		CODE		CODE		Ø inner	PRICE
	3/16" - 4,8 mm		1/4" - 6,4 mm		5/16" - 7,9 mm		3/8" - 9,5 mm		1/2" - 12,7 mm		
K0932-08		K0933-08		K0934-08		K0935-08		K0936-08		strong 8 oz - 224 g	\$4.20

Caution: these products contain natural rubber latex that may cause allergic reactions.

SAFETY RELEASE MODULES

For use with the suitable headgear or neck pads only. Safety clasps disengage the spring strap when tension force exceeds 2 Kg each side.

This system avoids dangerous slingshot rebound injuries of the facebow in case of disconnection from headgear tubes. Coil springs are incorporated into the modules to provide constant traction force. Available in three colour-coded attachment straps showing the force level.

Packs of 10

CODE		color	PRICE
M0810-01	light	transparent	\$71.10
M0810-02	medium	yellow	\$71.10
M0810-03	heavy	blue	\$71.10



STAINLESS STEEL SPRINGS


RELEASE PERMANENT FORCES

SAFETY RELEASE SPRING SYSTEM OF THE FORCE


NOTCH NUMBERED STRAPS FOR THERAPEUTICAL CONTROL

TOTAL SAFETY FOR THE PATIENT

CLASPS DISENGAGE THE SPRING STRAP WHEN TENSION FORCE EXCEEDS 2 KG EACH SIDE

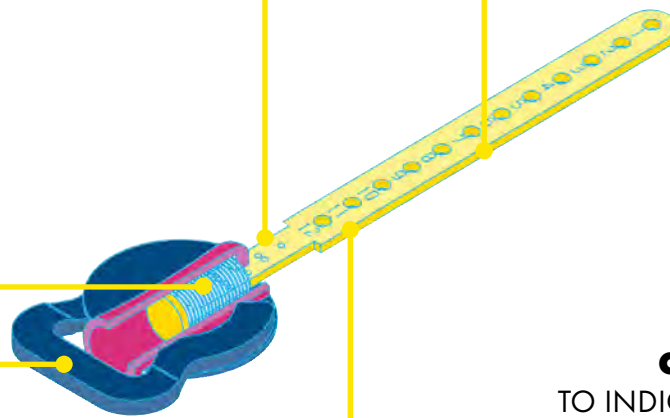

**DIFFERENTLY
COLORED STRAPS**
TO INDICATE THE DIFFERENT
SPRING FORCE OF MODULES

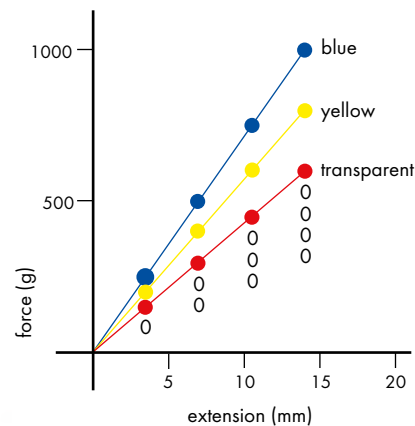

Traction force released by each module can be easily controlled: when 2 notch tension appears on the plastic strap under traction, 2/4 (two fourths) out of the total deliverable force of the spring is going to be released. Above diagram shows the pull forces released by each type of module.



EXTRAORAL FACEBOWS

Entirely made of biomedical spring stainless steel wire: the outer bow measures 1,80 mm in diameter and the single piece inner bow has a terminal portion of .045" (1,14 mm).

The laser welded bows joint, without the presence of solder material, is breakage resistant and avoids allergic reactions due to biometallism. Midline and number clearly laser marked.

Bright polished for maximum hygiene.

Outer bow is available in three lengths.

Standard size is available with pre-formed loops, bites or double bites, all with the option of hooks for elastics.

Caution: All extraoral facebows are intended for use with Leone M0810 safety modules.




STANDARD EXTRAORAL FACEBOWS

Thanks to the choice of making the inner stainless steel one-piece bow integral with the end portion of .045" (1,14 mm), the realization of loops and/or bayonet stops does not decrease its structural resistance, by ensuring the transmission of extra oral forces to the tubes on the bands.

Outer bow available in three lengths and with the option of elastic hooks.

Caution: All extraoral facebows are intended for use with Leone M0810 safety modules.

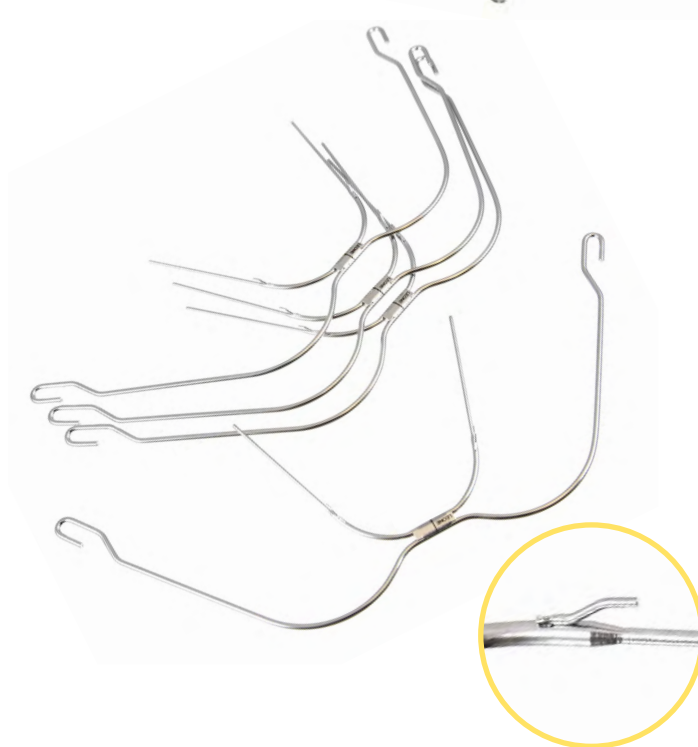


		CODE	pkg.	PRICE	CODE	pkg.	PRICE	outer arch diameter	inner arch diameter	
		standard facebow			standard facebow with hooks for elastics			mm	mm	inch
NEW	extralong outer arch	M4840-00	1	\$17.30				1,80	1,14	.045
		M4840-001	10	\$137.50				1,80	1,14	.045
	long outer arch	M4800-00	1	\$14.60	M4801-00	1	\$19.40	1,80	1,14	.045
		M4800-001	10	\$115.40	M4801-001	10	\$154.70	1,80	1,14	.045
	short outer arch	M4830-00	1	\$14.60	M4831-00	1	\$19.40	1,80	1,14	.045
		M4830-001	10	\$115.40	M4831-001	10	\$154.70	1,80	1,14	.045

EXTRAORAL FACEBOWS WITH PRE-FORMED LOOPS


Thanks to the choice of making the inner stainless-steel one-piece bow with the end portion of .045" (1,14mm) or .051" (1,3mm), the preformed loops are easily adaptable to the patient needs but also very breakage resistant. The .045" is available in five sizes while the .051" is available in three sizes.


Caution: all extraoral facebows are intended for use with Leone M0810 safety modules.

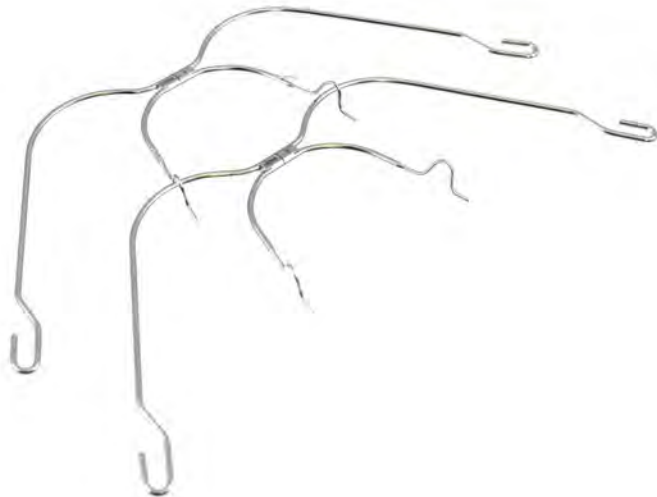

midline
clearly marked

WITH PREFORMED LOOPS				CODE	pkg.	PRICE	CODE	pkg.	PRICE	CODE	pkg.	PRICE	CODE	pkg.	PRICE
d = intermolar width from loop to loop in millimeter				d			long outer bow			long outer bow			short outer bow		
				83			M4805-01	1	\$17.30	M4805-011	10	\$137.00	M4835-01	1	\$17.30
				90			M4805-02	1	\$17.30	M4805-021	10	\$137.00	M4835-02	1	\$17.30
				97			M4805-03	1	\$17.30	M4805-031	10	\$137.00	M4835-03	1	\$17.30
				104			M4805-04	1	\$17.30	M4805-041	10	\$137.00	M4835-04	1	\$17.30
				111			M4805-05	1	\$17.30	M4805-051	10	\$137.00	M4835-05	1	\$17.30
													M4835-011	10	\$137.60
													M4835-021	10	\$137.60
													M4835-031	10	\$137.60
													M4835-041	10	\$137.60
													M4835-051	10	\$137.60
															size

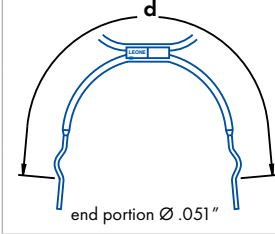WITH PREFORMED LOOPS				CODE	pkg.	PRICE	
d = intermolar width from loop to loop in millimeter				d			size
				90			2
				97			3
				104			4

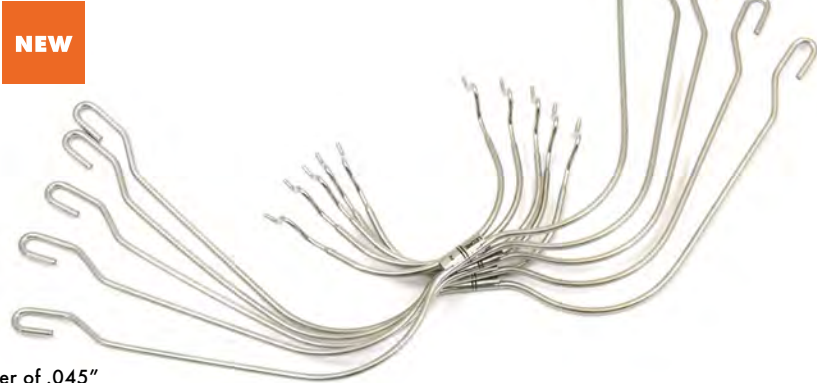

Caution: these are not intended to be used with the tubes with inner diameter of .045"

**M3309-00****FACEBOW SCREW STOP**

Made of stainless steel. For use with wire 1,14 mm (.045") in diameter.

Pack of 10

PRICE \$34.70

M3309-95**FACEBOW SCREW STOP KIT**

It includes:

50 stop screws and 1 screw-driver

PRICE \$185.60

**M0790-00****HEADGEAR WITH CHIN CAP**

For correction of Class III malocclusions with prognathism of the mandible. Headgear and elastic tractions are fully adaptable; the chin cap in natural leather assures an ideal fit to any face. Color: blue.

Pack of 1

PRICE \$52.00

REPLACEMENT PARTS**M0791-00**

Chin pad Pack of 1

PRICE \$25.20

M0791-01

Elastics Pack of 2

PRICE \$13.40

LIP BUMPER

Made of biomedical stainless steel .045" wire.
Non toxic techno-polymer labial coating.



M3050-00 without loops

PRICE \$93.50

M3051-00 with loops

PRICE \$94.00

**LIP BUMPER WITH LOOPS AND BAYONET STOPS**

CODE	pkg.	PRICE	CODE	pkg.	PRICE	size	mm
M3051-01	1	\$13.10	M3051-011	10	\$117.90	1	76
M3051-02	1	\$13.10	M3051-012	10	\$117.90	2	80
M3051-03	1	\$13.10	M3051-013	10	\$117.90	3	84
M3051-04	1	\$13.10	M3051-014	10	\$117.90	4	88
M3051-05	1	\$13.10	M3051-015	10	\$117.90	5	92



SPECIAL:
BUY 3
10-PACKS
GET ONE
FREE

CLASS II CORRECTOR

The non-compliance device for Class II correction.

The small sizes of the device allow for optimum patient comfort while the constant and light force delivered by the spring **MEMORIA®** 200 g, located inside the plunger, stimulates the mandibular advancement. In this **new version**, the **fluidity of movement and the strength of the telescopic mechanism are improved**. The packages include all the needful parts for the application of a bilateral correction device.



Available in **3 lengths** with both mesial and distal fitting to upper molar tube, for a total of 6 possible positions in the mouth. When used in the mesial position, it can be placed on a single direct-bonding tube without the need of band or round tube.




EXTRA-SHORT

MESIAL SUPPORT TO VESTIBULAR TUBE



DISTAL SUPPORT TO VESTIBULAR TUBE

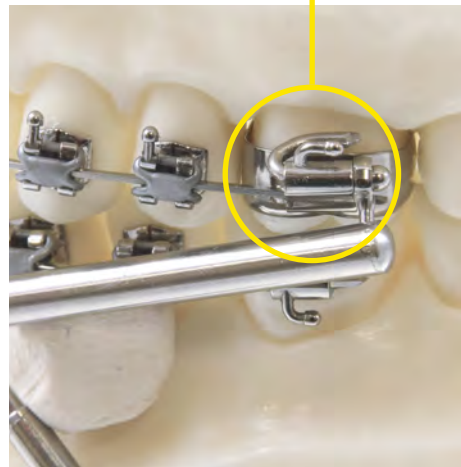

SHORT




LONG



1:1

To avoid breakage, the device should not be mounted at full compression: once in the mouth in habitual occlusion, a portion of the plunger coming out of the tube for at least **2 mm** should be visible.

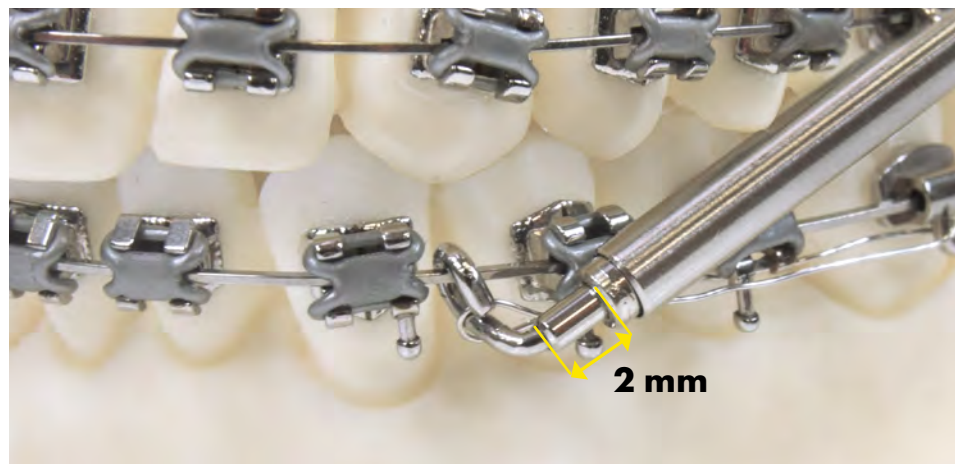

AVAILABLE IN 2 VERSIONS:

1 WITH ANTERIOR EYELET



2 WITH ANTERIOR HOOK



CODE	PRICE	CODE	PRICE					
WITH ANTERIOR EYELET		WITH ANTERIOR HOOK		Type	Length at full expansion (passive configuration)	Length at full compression (including 2 mm additional possible compensation)	Maximum safety stroke of spring	Maximum opening allowed to prevent disassembly
M2302-00	\$65.90	M2301-00	\$65.90	extra short	31 mm	25 mm	7,5 mm	48 mm*
M2302-01	\$65.90	M2301-01	\$65.90	short	35 mm	28 mm	10 mm	56 mm*
M2302-02	\$65.90	M2301-02	\$65.90	long	40 mm	33 mm	10 mm	62 mm*

* Distance between the 2 eyelets or hooks of anchorage

M2301-91

CLASS II CORRECTOR KIT WITH ANTERIOR HOOK

Pack of **3 correctors** with anterior hook (1 extra-short, 1 short, 1 long) and 6 measuring gauges

PRICE \$173.30

M2302-91

CLASS II CORRECTOR KIT WITH ANTERIOR EYELET

Pack of **3 correctors** with anterior eyelet (1 extra-short, 1 short, 1 long) and 6 measuring gauges

PRICE \$173.30



With the help of this measuring gauge, it will be simple and intuitive to determine which of the three available lengths will be the correct one.



P1904-00 **LIGHT WIRE PLIERS**

Longer and more tapered beaks than Angle pliers. Excellent for bending small diameter loops. For wires up to .018" (0,46 mm).

PRICE \$111.80



P1154-00 **OPTICAL STYLE** **PLIER, MEDIUM TIP**

For wires up to .028" (0,71 mm).

PRICE \$121.00



P1155-00 **OPTICAL STYLE** **PLIER, LARGE TIP**

For wires up to .032" (0,80 mm).

PRICE \$100.30



P1905-00 **RECTANGULAR ARCH** **PLIERS TWEED STYLE** **SHORT JAW**

For use with square or rectangular wires. The blades are designed to be parallel at .020" (0,51 mm) pliers opening. Working edges are carefully bevelled with a diamond hone to prevent wire scoring.

PRICE \$129.70



P1906-00 **THREE JAW PLIERS** **ADERER STYLE MEDIUM**

For wires up to .036" (0,91 mm).

PRICE \$111.80



P1075-01 **THREE JAW PLIERS** **ADERER STYLE THIN**

For round wires up to .020" (0,51 mm) and rectangular up to .021"x.025" (0,53x0,64 mm).

PRICE \$217.60



P1900-00
BIRD BEAK PLIERS
ANGLE STYLE

Necessary for working round wires up to .032" (0,80 mm) in diameter. The beaks are .040" (1,02 mm) at the tip. A 9° angle is ground on the sides of pyramid beak to allow for wire springback. Beaks are parallel at .020" (0,51 mm) pliers opening.

PRICE \$181.90



P1901-00
LOOP FORMING PLIERS
TWEED STYLE

A plier designed to make loops in round or rectangular wires up to .022"x.028" (0,56x0,71 mm). The cylindrical beaks have two sections: .043" (1,09 mm) and .063" (1,60 mm).

PRICE \$197.90



P1903-00
YOUNG STYLE PLIERS

For wires up to .040" (1,02 mm).

PRICE \$126.00



P1907-00
HOW STYLE PLIERS
STRAIGHT

The 3/32" (2,4 mm) gripping pads are serrated and shaped so that their center line coincides with the axis of rotation of the pliers. Beaks are rounded to be completely lip safe.

PRICE \$111.80



P1908-00
HOW STYLE PLIERS
STRAIGHT

Exactly the same as P1907-00. With 1/8" (3,2 mm) gripping pads.

PRICE \$111.80



P1909-00
HOW STYLE PLIERS
OFFSET

Exactly the same as P1908-00 with tips incorporating a 45° offset angle.

PRICE \$129.70



P1910-00 WEINGART STYLE PLIERS

Classical design, the tips are carefully rounded to make this pliers completely mouth safe. Useful for removal of ceramics.

PRICE \$129.70



P1915-00 DISTAL END CUTTER WITH SAFETY HOLD

The small dimensions do not interfere with quality. It cuts archwires and safely holds the loose distal end.

Caution: maximum cutting capacity:

- rectangular wire up to .022"x.028" (0,56x0,71 mm);
- round wire up to .020" (0,51 mm);
- twist wire up to .016" (0,41 mm).

PRICE \$226.50



P1911-00 STOP AND SURGICAL HOOK CRIMPING PLIERS

A plier designed to crimp surgical hooks in the archwire. It can also be used to fix round or rectangular stops.

PRICE \$143.30




P1914-00 MICRO DISTAL END CUTTER WITH SAFETY HOLD LONG HANDLES

The retracted position of 10° angled hand grip affords a wider view in the oral cavity.

Caution: maximum cutting capacity:

- rectangular wire up to .022"x.028" (0,56x0,71 mm);
- round wire up to .020" (0,51 mm);
- twist wire up to .016" (0,41 mm).

PRICE \$253.60

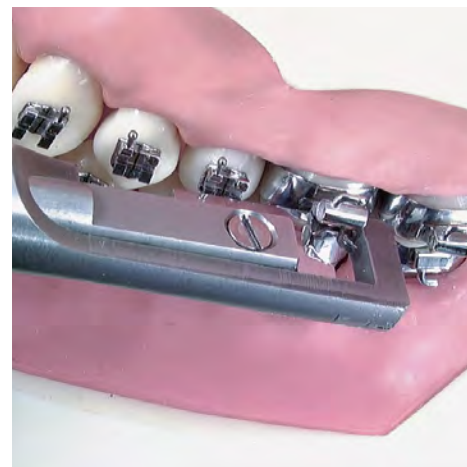

P1090-00 CAP REMOVER FOR CONVERTIBLE TUBES

It makes the conversion of rectangular tube easy when banding of the second molar is required. Autoclavable.

PRICE \$222.10

P1090-01 REPLACEABLE TIPS P1090-00

Pack of 2 **PRICE \$20.50**



P1530-00
LABORATORY PLIERS
ADAMS STYLE

Large tips with tungsten carbide inserts. Necessary for bending wires in the laboratory. A round cone tip of .052" (1,3mm) and a pyramid beak of .056" (1,4mm). For wires up to .060" (1,5mm)

Excellent for Crozat technique.

PRICE \$195.30



P1034-00
THREE JAW PLIERS
ROUND

Very sturdy, does not score wires. Suitable for bending the arms of expanders. For wires up to 1,5 mm. (.059")

PRICE \$175.10



P1099-00
HEAVY DUTY CUTTER

For hard wires up to .060" (1,5mm)

PRICE \$147.00



P1162-00
UNIVERSAL PLIERS

Extremely tapered design. For wires up to .036" (0,9mm)

PRICE \$44.90



P1163-00
UNIVERSAL PLIERS

Classic design. For wires up to .036" (0,9mm)

PRICE \$44.40



P1917-00
HARD WIRE CUTTER

For cutting orthodontic wires for fixed technique any size.

Caution: maximum cutting capacity:
 - rectangular wire up to .022"x.028" (0,56x0,71 mm);
 - round wire up to .020" (0,51 mm).

PRICE \$195.30



P1919-00
MINIATURE
LIGATURE CUTTER

For soft wires.

Caution: maximum cutting capacity:
ligature wire up to .014" (0,36 mm).

PRICE \$228.60




P1923-00
POSTERIOR BAND
REMOVING PLIERS
SHORT

PRICE \$159.10

Necessary for removal of posterior bands. A plastic padded tip and a sharp removing beak allow, with a slight pressure, high precision band removal.

P1453-00
REPLACEMENT PAD
3/16" (4,8 MM)
FOR P1923-00

PRICE \$15.80

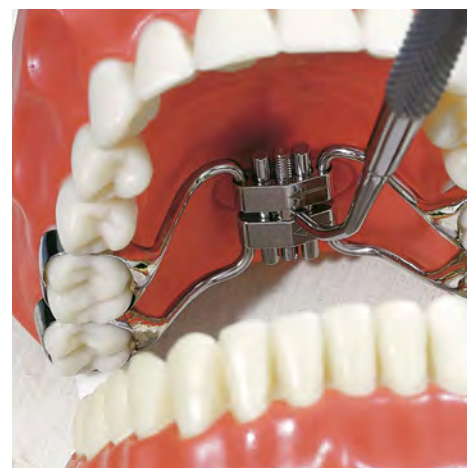

P1621-00
INSTRUMENT FOR
ENDORAL ACTIVATION

Suggested for use with multifunctional orthodontic screws.

PRICE \$18.90



P1509-00
STAINLESS STEEL UNIVERSAL
BRACKET ALTIMETER

A small device for measurement of the right height of bracket slots.

PRICE \$81.60



P1026-00
D.B. BRACKET TWEEZERS

It holds any type of bracket, from a wide twin to a single one.

Self-locking of the bracket and release by pressure.

Several tweezers can be pre-loaded with brackets to bond them one after the other.

PRICE \$40.40

P1098-00**UTILITY INSTRUMENT****PRICE \$47.00****P1139-00****LIGATURE TUCKER**

Double ended for placement or removal of elastics. Necessary for tucking metal ligatures. **PRICE \$39.10**

**P1001-01****BAND DRIVER OFFSET TIP**

Fully made of stainless steel.
For placement and adjustment of bands. Activation is ease in all positions.
Seating pressure range 80 up to 1000 g.

PRICE \$236.30**P1002-01****REPLACEMENT TIP
FOR P1001-01****PRICE \$35.70****P1000-01****BAND DRIVER STRAIGHT TIP**

Exactly the same as P1001-01 but with straight tip.

PRICE \$236.30**P1003-01****REPLACEMENT TIP
FOR P1000-01****PRICE \$35.70****P1505-01****TWISTER****FOR PREFORMED LIGATURE TIES
K4430-25 AND K4480-30**

Fully made of stainless steel. Push the back button to insert the preformed ligature tie. It makes position and final legation easy.

PRICE \$153.30



P1167-00
TURRET
BLUE CODE

For forming rectangular arch wires
.016" up to .0215".

PRICE \$199.20



P1168-00
TURRET
GOLD CODE

For forming rectangular arch wires
with torque adjustments: 0° - 7° - 10° -
13° - 16°. For use only with wire: .018".

PRICE \$199.20



LIGATURE TYING PLIER
MATHIEU STYLE

P1160-00 Tips of 2,5 mm
PRICE \$71.10

P1161-00 Tips of 1,5 mm
PRICE \$69.80






P1089-00
MOSQUITO
FORCEPS OFFSET

PRICE \$8.10

P1088-00
MOSQUITO
FORCEPS STRAIGHT

PRICE \$8.10

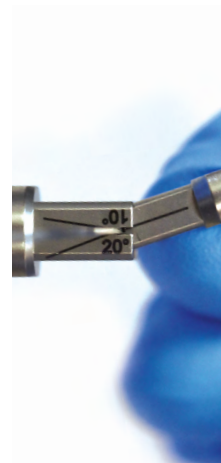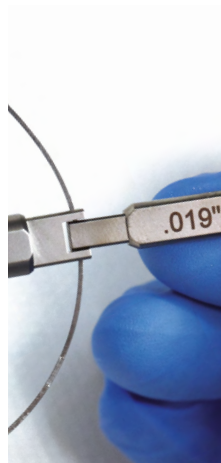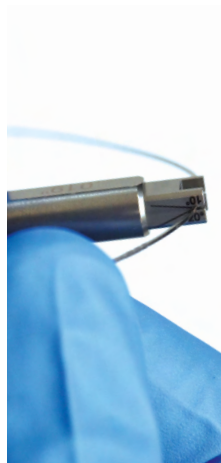

P1635-00
TORQUING INSTRUMENT
FOR WIRE .019"X.025" AND .021"X.025"/.027"

A device consisting of two components to easily apply torque to a portion of arch wire in proximity of a specific bracket. The instrument with the U end is used to mark out and maintain the selected portion of the arch wire. The additional instrument is used to firmly keep the rectangular wire and to apply the desired torque degrees. Thanks to the presence of a laser marked protractor, it is possible to check the degrees of torque given. One instrument end is for arch wires .019"x.025", the other for .021"x.025"/.027"

PRICE \$124.20





R6000-00 LEONEAMERICA COLOR GUIDE

Colors illustrated are only indicative. Color may slightly change compared to those obtained with **LEONEAMERICA®** acrylic resin.

PRICE \$13.70



R3712-00 QUADRANGULAR WAX STICKS

Wax in 25 cm long sticks, 3x3 mm cross section. The special adherence and high workability make it useful to border the impression trays to fit the edges of the functional impression.

Pack of 110 sticks

PRICE \$60.40



R3710-00 UTILITY WAX

Highly soft wax in strips which bends at room temperature without breaking. Thanks to its good adherence and high workability, it serves to border the impressions to obtain the base of the mould.

Pack of 24 strips

PRICE \$60.40



R3711-00 MODELLING WAX

Wax in sheets. It moulds, both cold and hot, without breaking.

Pack of 500 g

PRICE \$47.30



R3730-00 BITE WAX

Necessary for taking accurate occlusal impressions. A special aluminium foil lessens distortions, 5 mm thick.


Pack of 24

PRICE \$84.80


BIOCOMPATIBLE THERMOFORMING MATERIAL


Discs and square plates made out of PET G (Polyethylene Terephthalate modified with Glycol), a material ideal for thermoforming to manufacture dental aligners and post-treatment aesthetic retainers.


The material conforms to ISO 10993-1, and once formed, it features excellent optical transparency, with highly efficient biomechanical properties.

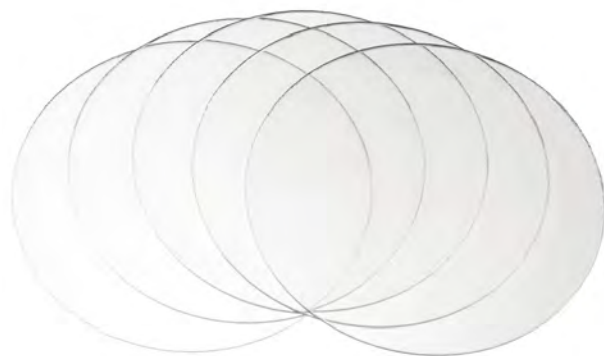

**CALL TO
RECEIVE A FREE
SINGLE DISC OR PLATE
SAMPLE**

THERMOFORMABLE DISCS Ø 120 mm	CODE	thickness	pkg.	PRICE
	R7001-05	0,5 mm	25 pcs	\$54.30
	R7001-08	0,8 mm	25 pcs	\$54.30
	R7001-10	1,0 mm	25 pcs	\$54.30
	R7001-15	1,5 mm	15 pcs	\$35.80
	R7001-20	2,0 mm	10 pcs	\$22.80

THERMOFORMABLE DISCS Ø 125 mm	CODE	thickness	pkg.	PRICE
	R7002-05	0,5 mm	25 pcs	\$54.30
	R7002-08	0,8 mm	25 pcs	\$54.30
	R7002-10	1,0 mm	25 pcs	\$54.30
	R7002-15	1,5 mm	15 pcs	\$35.80
	R7002-20	2,0 mm	10 pcs	\$22.80

THERMOFORMABLE PLATES 125 x 125 mm	CODE	thickness	pkg.	PRICE
	R7013-05	0,5 mm	25 pcs	\$54.30
	R7013-08	0,8 mm	25 pcs	\$54.30
	R7013-10	1,0 mm	25 pcs	\$54.30
	R7013-15	1,5 mm	15 pcs	\$35.80
	R7013-20	2,0 mm	10 pcs	\$22.80

**SPECIAL:
BUY 4 GET 1 FREE OF EQUAL
OR LESSER VALUE
MIX & MATCH**


BIOCOMPATIBLE THERMOFORMING MATERIAL HARD/SOFT


NEW


Discs and square plates featuring one soft and one hard side. The soft side is made of EVA (Ethylene Vinyl Acetate) and the hard side is made of PET-G (Polyethylene Terephthalate modified with Glycol). Acrylic resin is perfectly combined with the hard side of the thermo-formable material, making it suitable for the fabrication of splints for people suffering from bruxism or as a retentive base in the construction of complex acrylic appliances, such as anti-snoring devices. Once thermoformed, the material will get a high breakage resistance.



**CALL TO
RECEIVE A FREE
SINGLE DISC OR PLATE
SAMPLE**

HARD/SOFT DISCS Ø 120 mm	CODE	thickness	pkg.	PRICE
	R7101-02	2 mm	10 pcs	\$119.00
	R7101-03	3 mm	10 pcs	\$149.00

HARD/SOFT DISCS Ø 125 mm	CODE	thickness	pkg.	PRICE
	R7102-02	2 mm	10 pcs	\$119.00
	R7102-03	3 mm	10 pcs	\$149.00

HARD/SOFT PLATES 125x125 mm	CODE	thickness	pkg.	PRICE
	R7113-02	2 mm	10 pcs	\$119.00
	R7113-03	3 mm	10 pcs	\$149.00



SPECIAL:
**BUY 4 GET 1 FREE OF EQUAL
OR LESSER VALUE
MIX & MATCH**

CALIBRA® BANDS
UPPER BICUSPIDS
CALIBRA® BANDS
LOWER BICUSPIDS

5 4 | 4 5

as indication only

5 4 | 4 5

LEONE Calibra - E5101/E5100-00										LEONE Calibra - E5501/E5500-00									
G.A.C. Dentsply Sirona Snap-Fit*	AMERICAN ORTH. Mr. Bicuspid*	ORMCO Bicuspid Bands*	DENTAURUM Standard*	UNITEK Contoured*	UNITEK General purpose*	UNITEK Proportioned*	FORESTADENT*	R.M.O. Truform*	LANCER Universal*	G.A.C. Dentsply Sirona Snap-Fit*	AMERICAN ORTH. Mr. Bicuspid*	ORMCO Bicuspid Bands*	DENTAURUM Standard*	UNITEK Contoured*	UNITEK General purpose*	UNITEK Proportioned*	FORESTADENT*	R.M.O. Truform*	LANCER Universal*
1		1	1		00		1			1			1				1		
2	1										1		2				2		1
3	2	2	2				2	G	1	2		1	3		0	1	3		
4	3	2	3		0		3			3	2		4		1+	4	4	F	2
5	4	3	4		0+	1	4	G+	2	4		3	5	1	2	5	5		
6			5		1	2	5			5	5	4	6		2+	6	6	F+	3
7	1	5	4	6	1+		6	H	3	6	6	6	7		3+	8	7		4
8		7	6	7	2+	4	7	H+		7	8								
9	2	8		8	2		8		5	8									
10	3	9	7	9		3	6	9	6	9	1	9	7		5	10	8	G	6
11	4		9	10		4		10	7			8	9			11	9	G+	7
12		11	10	11		4+	8	11	I+		2	10	9	10	5	5+	12	10	H
13	5	12	11	12	4		9	12	8	10	3	11	10	11		6+	11		8
14	6	14	12	13		5	10	13	J	11			11	12	7	7	14	11	H+
15	7	15	13	14		6	12	14	J+	12	4	13		13		8	16	13	10
16	8	17	14	15	6	6+		15	11	13		14	12	14	9	8+	17	14	I
17				16		7	13	16	K	14									
18	9	19	15	17	7		14	17		15	5	15		15		9+		15	I+
19	10	20	16	18	8	7+	15	18		16	6	16	14	16	11	10	19	16	13
20			17	19		8	16	19	K+	17			15	17	12	11	21	17	J
21		21								18	7	17	17	18		12		18	15
22		22								19	9	19		19	14	12+	23	19	K
23	12	24	19	20	10	9+		20	L	20									
24	13	25	20	21		10	19	21		21	10	20	19	20	15	13+	25	20	
25	14		21	22			20	22	L+	22	11	21		21		14+	27	21	K+
	15	27	22	23	11	10+	21	23	M	23	12	23	21	22	17	15		22	18
26		28	23	24		11	22	24		24		24		23	18	15+	28	23	L
27	16		24	25	12	11+	23	25		25	13		22	24		16	29	24	L+
28		29								26		26							
29		30	25	26		12+	24	26	M+	27		28		25		17	31	25	M
30	17		26	27	13	13	25	27		28	14	29	25	26		18	33	26	22
31		31			14					29									
32	18	32	28	28		13+	27	28	N	30		30	26	27		19	34	27	M+
			29	29		14	28	29	N+	31	16		27	28	22	19+	35	28	
	19		31	30		15	29	30		32	17	31	28	29		20	36	29	24
			32	31			30	31				32	29	30				30	25
	20			32				32					30	31			38	31	26
													32					32	27

* it is not a trademark of Leone S.p.A. and Leone America Dental Products inc.

CALIBRA® UPPER 1 st MOLAR														CALIBRA® LOWER 1 st MOLAR																					
6		6		as indication only																												6		6	
LEONE Calibra - E1200/E1300-00														LEONE Calibra - E1600/E1700-00																					
G.A.C. Dentsply Sirona Snap-Fit*														G.A.C. Dentsply Sirona Snap-Fit*																					
AMERICAN ORTH. MRI Contoured*														AMERICAN ORTH. MRI Contoured*																					
ORMCO Standard High-Retention*														ORMCO Standard High-Retention*																					
ORMCO Ultima*														ORMCO Ultima*																					
ORMCO Washbon*														ORMCO Washbon*																					
DENTAURUM Standard*														DENTAURUM Standard*																					
DENTAURUM Dentaform*														DENTAURUM Dentaform*																					
UNITEK Narrow Contoured*														UNITEK Narrow Contoured*																					
UNITEK General Purpose*														UNITEK General Purpose*																					
FORESTADENT W Fit Form*														FORESTADENT W Fit Form*																					
HENRY SCHEIN Epic*														HENRY SCHEIN Epic*																					
R.M.O. Truform*														R.M.O. Truform*																					
LANCER Original*														LANCER Original*																					
1														1																					
2														2																					
3														3																					
4														4																					
5														5																					
6														6																					
7														7																					
8														8																					
9														9																					
10														10																					
11														11																					
12														12																					
13														13																					
14														14																					
15														15																					
16														16																					
17														17																					
18														18																					
19														19																					
20														20																					
21														21																					
22														22																					
23														23																					
24														24																					
25														25																					
26														26																					
27														27																					
28														28																					
29														29																					
30														30																					
31														31																					
32														32																					
33														33																					
34														34																					
35														35																					
36														36																					

* it is not a trademark of Leone S.p.A. and Leone America Dental Products inc.

CONVERSION CHARTS

inch	mm	inch	mm	inch	mm	inch	mm	inch	mm
.001	0,03	.017	0,43	.033	0,84	.049	1,24	.065	1,65
.002	0,05	.018	0,46	.034	0,86	.050	1,27	.066	1,68
.003	0,08	.019	0,48	.035	0,89	.051	1,30	.067	1,70
.004	0,10	.020	0,51	.036	0,91	.052	1,32	.068	1,73
.005	0,13	.021	0,53	.037	0,94	.053	1,35	.069	1,75
.006	0,15	.022	0,56	.038	0,96	.054	1,37	.070	1,78
.007	0,18	.023	0,58	.039	0,99	.055	1,40	.071	1,80
.008	0,20	.024	0,61	.040	1,02	.056	1,42	.072	1,83
.009	0,23	.025	0,64	.041	1,04	.057	1,45	.073	1,85
.010	0,25	.026	0,66	.042	1,07	.058	1,47	.074	1,88
.011	0,28	.027	0,69	.043	1,09	.059	1,50	.075	1,90
.012	0,30	.028	0,71	.044	1,12	.060	1,52	.076	1,93
.013	0,33	.029	0,74	.045	1,14	.061	1,55	.077	1,96
.014	0,36	.030	0,76	.046	1,17	.062	1,57	.078	1,98
.015	0,38	.031	0,79	.047	1,19	.063	1,60	.079	2,01
.016	0,41	.032	0,81	.048	1,22	.064	1,63	.080	2,03

inch	mm
1/64"	= .0156
1/32"	= .0312
1/16"	= .0625
1/8"	= .125
1/4"	= .250
3/8"	= .375
1/2"	= .5
5/8"	= .625
3/4"	= .75
7/8"	= .875
1"	= 1.0

1 mg	=	0,0154 grains
1 g	=	1000 mg = 0,0353 oz
1 Kg	=	1000 g = 2.2046 lb
1 ounce	=	28,35 g
1 pound	=	16 ounce = 0,4536 Kg

1 inch = 25,40 mm
1 mm = 0,0394 inch

1 N/mmq = 145.14 p.s.i.
1 p.s.i. = 0,00689 N/mmq

FDI SYSTEM	18	17	16	15	14	13	12	11	21	22	23	24	25	26	27	28
	48	47	46	45	44	43	42	41	31	32	33	34	35	36	37	38

CATALOGUE SYMBOLS

On the various tables of catalogue we used symbols to better indicate the characteristics of our products. Below are the symbols used:	dimensions		angulation	
	diameter		rotation of tube	
one turn of screw	inner diameter		upper right	
expansion	torque of bracket		upper left	
length	torque of tube		lower right	
height	slot		lower left	

Safety data sheets of our products can be downloaded from our web site: www.leone.it

PRODUCT LABEL SYMBOLS

The label on the package of any medical device set on the market will show the symbols in compliance with the harmonized standards. The symbols marked with a single (*) are based on the ISO 21531, ISO 15223-1 and on the 93/42EEC Directive. The symbols marked with double (**) have been instead adapted by Leone.

manufacturer's trade name and address		manufacturer's catalogue code		code number and product description in different languages	
CE mark (made in compliance with 93/42EEC Directive on class IIA or IIB medical devices)		expiry date, if the product is perishable (year/month/day)		storage temperature	
lot number * (indicated by LOT mark)		for professional use only		for single use only	
keep dry		this product contains Nickel-Chromium: possible allergic reactions		keep away from sunlight	
CE mark (made in compliance with 93/42EEC Directive on Class I medical devices)		read the documents enclosed		gamma-ray sterilized	
titanium		surgical steel		this product contains Chromium: possible allergic reactions	
autoclavable at temperature indicated		with content or presence of natural rubber latex		non-sterile	
do not use if package is damaged		read the instructions for use enclosed		bar code or Data Matrix	

* Unless otherwise indicated, the LOT number explains the date of manufacture with the key YYMMDDNN (YY year, MM month, DD day, NN Leone progressive internal number of the lot). For example: 19021902 indicates the lot no. 2 of February 19, 2019.